

# DENTAL WORLD

## CATALOGUE 3°



OUR PRODUCTS YOUR SUCCESS!!



## CONSERVATIVE

FILL DENT	LIGHT-CURED RESTORATIVE NANO COMPOSITE RESIN	A 1
FILL DENT FLOW	LIGHT-CURED RADIOPAQUE FLOWABLE COMPOSITE RESIN	A2
ORTHO FLOW	LIGHT CURING ORTHODONTIC FLOWABLE COMPOSITE	A3
ULTRA ONE BOND	ONE COMPONENT SELF-ETCHING LIGHT-CURED ADHESIVE	A4
BASIC BOND	SINGLE COMPONENT BONDING AGENT	A4
FILL DENT SEAL	LIGHT CURED PIT & FISSURE SEALANT	A5
LIGHT TEMP	LIGHT-CURED TEMPORARY FILLING MATERIAL	A5
FILL ETCH	ETCHING AGENT WITHOUT SILICA	A6
OVER ETCH	NEW AND TOTAL ETCHING	A6
DUAL BUILD CORE	DUAL-CURED CORE BUILD-UP MATERIAL	A7
SELF BILD CORE	SELF-CURED CORE BUILD-UP MATERIAL	A7
ONE BUILD CEM	SELF-ADHESIVE UNIVERSAL RESIN CEMENT	A8
CANAL PREP	CREAMY EDTA ROTH CANAL CLEANER	A8
CANAL PASTE	CALCIUM HYDROXIDE PASTE	A8

## IMPRESSION

FLUENT ELASTIC	IMPRESSION MATERIAL PUTTY	B 1
FLUENT SOFT NORMAL	IMPRESSION MATERIAL PUTTY	B 1
FLUENT SOFT FAST	IMPRESSION MATERIAL PUTTY	B2
FLUENT HARD	IMPRESSION MATERIAL PUTTY	B2
FLUENTB LIGHT	IMPRESSION MATERIAL CARTRIDGE	B3
FLUENT REGULAR	IMPRESSION MATERIAL CARTRIDGE	B3
FLUENT OCCLUCAD	IMPRESSION MATERIAL CARTRIDGE	B4
FLUENT CHECK	IMPRESSION MATERIAL CARTRIDGE	B4
FLUENT CLEAR TIXO	IMPRESSION MATERIAL CARTRIDGE	B5
FLUENT IMPLA GUM	SEMI-RIGID TYPE	B5
DABLIU MIX	MIXING GUN	B5
DABLIU SIL	IMPRESSION MATERIAL	B6
DABLIU GUM	IMPRESSION MATERIAL	B6
DABLIU CATALIST	PASTE / LIQUID CATALYST FOR SILICONE RESINS	B6
DABLIU SIL LIGHT	IMPRESSION MATERIAL	B6
DABLIU SIL REGULAR	IMPRESSION MATERIAL	B6
DABLIU ALGIN	HIGH PRECISION CHROMATIC ALGINATE	B7
DABLIU ALGIN GREEN	HIGH PRECISION NON-CHROMATIC ALGINATE	B7
ALGINBOY	ORTHODONTIC ALGINATE	B7

## ENDODONTICS

ENDO PERFECT PAPER	ENDODONTIC PAPER POINTS	C 1
ENDO PERFECT PAPER	ENDODONTIC PAPER POINTS TAPER SIZE	C2
ENDO PERFECT GUTTA PERCHA	ENDODONTIC GUTTA PERCHA POINTS	C2
DABLIU FILES NI-TI	NI-TI CANAL INSTRUMENTS	C3
DABLIU FILES STEEL	STAINLESS STEEL CANAL INSTRUMENTS	C4
PASTE CARRIER	PASTE CARRIER	C5
BARBED BROACH	BARBED BROACH	C5
ENDO ACCESSORIES	ENDO CLEAN, ENDO PLASTIC BOX, ENDO BLOC GAUGE, SILICON STOPS	C5
DABLIU ENDO POST	ENDOCANAL POSTS	C6
DABLIU CONIX POST	ENDOCANAL POSTS	C7
DABLIU CROWN POST	ENDOCANAL POSTS	C7
SHORT POST R&R MICRO	ENDOCANAL POSTS	C8

## BLEACHING

DABLIU SMILE	BLEACHING GEL 38% (PEROXIDE HYDROGEN)	D 1
DABLIU DAM	LIGHT CURING LIQUID DAM	D 1
BRILLIANCE PEN	BLEACHING (CARBAMIDE PEROXIDE)	D2
BRILLIANCE	BLEACHING (CARBAMIDE PEROXIDE)	D3
THERMOFORMING	MACHINE	D3

## HANDPIECES

HANDPIECE TS-46 I	E-GENERATOR INTEGRATED LED HANDPIECE	E 1
HANDPIECE TS-26 I	HANDPIECE WITH AND WITHOUT FIBER OPTIC	E2


HANDPIECE TS-22 I	DISPOSABLE HANDPIECE	E2
HANDPIECE TS-42 I	RAINBOW HANDPIECE SERIES	E3
LOW SPEED TS SERIES	LOW SPEED AIR MOTOR	E3
LOW SPEED TS SERIES	LOW SPEED CONTRA ANGLE	E4
DISINFECTION		
SPORIDOX	DISINFECTANT	F1
SPOREXIN PLUS DETERSONIC	DISINFECTION BY ULTRASOUND	F1
PROTEOZIM PLUS 400	DISINFECTION BY ULTRASOUND	F1
PROTEAZONE OD	PRE-STERILIZATION	F2
ADASPOR READY	STERILIZATION	F2
SPOREXIN PLUS OD	STERILIZING AND HIGH LEVEL DISINFECTANT SOLUTION	F2
SPOREXIN PLUS DS	SPRAY FOR SURFACES	F3
SPOREXIN PLUS DS TOWELS	TOWELS FOR SURFACES	F3
ISACLEAN	SPRAY FOR SURFACES	F3
SPOREXIN PLUS VACUUM	DISINFECTANT FOR SOLUTION DEVICES	F4
GERMICIDIN	DISINFECTANT, DETERGENT, CONCENTRATED SOLUTION	F4
BLUSTERIL ALCOHOLIC	DISINFECTANT, CONSERVANT	F4
BENODENT	MOUTH WASH	F5
SPOREXIN	HAND SOAP	F5
DISINFECTION TRAY	DISINFECTION TRAY 1 LT. - 3LT.	F5
HAND SANITIZER GEL	HAND SANITIZER GEL	F6
VARIOUS		
ORTHO BOX & ORTHO BITE	ORTHO BOX & ORTHO BITE	G1
PROSTHESES CASE	PROSTHESES CASE	G1
DABLIU BRUSH	DABLIU BRUSH	G2
DABLIU ULTRA BRUSH	DISPOSABLE BRUSH APPLICATORS	G2
DABLIU ENDO BRUSH	ENDODONTIC MICRO BRUSH	G2
MICRO APPLICATORS	MICRO APPLICATORS FOR BRUSH HANDLE	G3
VISORY GALAXY	VISORY GALAXY	G3
PROTECTIVE GLASSES	PROTECTIVE GLASSES	G3
PROTECT PLUS	DENTAL BIBS WITH POCKET	G4
DABLIU TIPS	APPLICATORS NEEDLES	G4
DISPENSER FOR COTTON	DISPENSER FOR COTTON	G5
CRISTAL TOOTH	CRISTAL TOOTH	G5
APPLYTRAY	FLUORIDE GEL TRAYS	G5
DABLIU MIXING CUPS	DABLIU MIXING CUPS	G6
DABLIU CARPET	BACTERIOSTATIC CARPET	G6
DABLIU PROTECTOR	BARRIER FILM	G6
ANTISTATIS HEALTH CLOGS	COLORLED CLOGS	G7
SLEEVES	X-RAY BARRIER ENVELOPES	G8
SLEEVES	X-RAY	G8
SLEEVES	SUCTION TUBE	G8
SLEEVES	CHAIR & HEADREST	G8
SLEEVES	X-RAY SENSOR & MOUTH MIRROR	G9
SLEEVES	X-RAY SENSOR	G9
SLEEVES	INTRAORAL CAMERA	G9
SLEEVES	T-STYLE LIGHT HANDLE	G9
SLEEVES	PEN TYPE CURING LIGHT	G10
SLEEVES	CURING LIGHT	G10
SLEEVES	SYRINGES	G10
SLEEVES	HANDPIECE	G10
LABORATORY		
DABLIU LAB	POLYCONDENSATION VINYL POLYSILOXANE	H1
DABLIU LAB A+B	ADDITION VINYL POLYSILOXANE HIGH PRECISION IMPRESSION	H1
FLUENT CLEAR FLUID	IMPRESSION REALIZATION FOR MASKS	H2
FLUENT GUM & SCAN	IMPRESSION GINGIVAL REPRODUCTION	H2
DUPLEX LAB 1 4-25	SILICONE	H3
DUPLEX ANGLE 1 2-22-27	SILICONE	H3
DABLIU WAX	RED WAX	H4
DABLIU PUMICE	NATURAL POWDER PUMICE	H4
DABLIU STONE IV-EH	GYPSUM TYPE IV GOLDEN BROWN	H4
DABLIU STONE MODEL IV	GYPSUM TYPE IV FOR MODEL	H4
DABLIU STONE ORTHO III	GYPSUM TYPE III ORTHODONTIC WHITE	H5
DABLIU STONE ORTHO IV	GYPSUM TYPE IV ORTHODONTIC WHITE	H5
DABLIU STONE MODEL III	GYPSUM TYPE III FOR MODELS	H5
DABLIU STONE MODEL II	GYPSUM TYPE II FOR MODELS	H5





GINGIVAL RETRACTION & PERFECT BLEEDING CONTROL

**Plus**  
*GingiMaster*

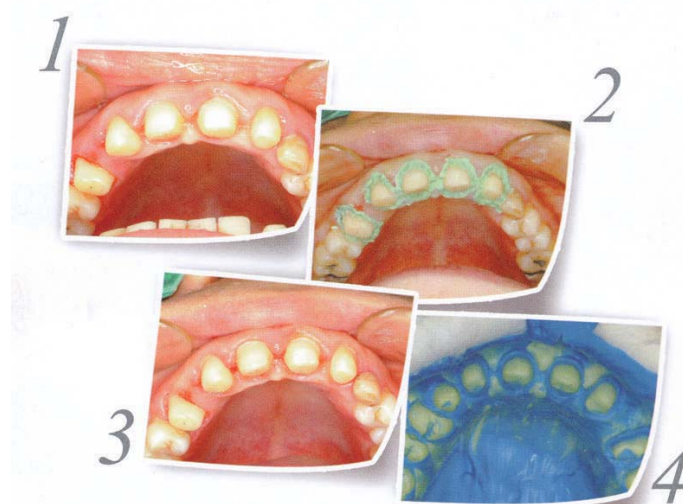

## EFFECTIVE

- Impression taking
- For restoring class II or V cavities
- For sealing prostheses
- For cementations

## SIMPLE & RAPID

- GingiMaster is injected into the sulcus using the applicator
- Only 15 to 20 sec required for product application
- Sulcular opening is obtained in 30 sec to 2 min depending on gingival texture
- The end-result is a dry, wide-open sulcus with access to the margins.

## SAFE

- There is no danger of gingival recession or bone loss related damage of the Epithelial Attachment
- GingiMaster Plus exerts a moderate calculated pressure on the marginal gingival, slightly higher than the pressure the gum exerts on the tooth.

## PAINLESS

- No anaesthesia is required with GingiMaster Plus
- Complete comfort for the patient

**EXCLUSIVE**  
DEALER FOR EUROPE



# FILL DENT

LIGHT-CURED RESTORATIVE NANO COMPOSITE RESIN



**Fill Dent** is a light-cured restorative nano-hybrid composite resin indicated for both Posterior and Anterior restorations. It is designed for minimizing polymerization shrinkage by inorganic fillers high loading, thus **Fill Dent** has excellent marginal sealing and surface texture. By adopting optimal polymerization initiator system, **Fill Dent** has high strength and abrasion resistance. It contains radiopaque filler; it also has an optimal viscosity for restoration and shows an excellent workability.

## Capsule Dispenser


- Ergonomic;
- Autoclavable.




007721 | Capsule Dispenser

## Indications:


- Restorations for Posterior and Anterior
- Core Build Up
- Splinting

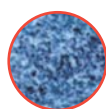

Fill Dent

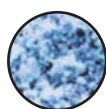

Company D

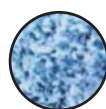

Company K

## Packages:

### Kit

- 4g x 5syringes A1, A2, A3, A3.5, B2)
- Bonding agent 5ml x 1bottle
- Etchant 3ml x 1syringe
- Brushes, Brush holder Disposable tips(blue), accessories

## Refill Syringe




4gx 1syringe

- Shade: A1, A2, A3, A3.5, A4, B1, B2, B3, C2, C3, D3, A2O, A3O, EC Enamel Cervical, T Transparent, EI Enamel Incisal

## Refill Capsule

0.25g x 20ea

- Shade: A1, A2, A3, A3.5, A4, B1, B2, C3, D3, A2O, A3O, EC Enamel Cervical, T Transparent, EI Enamel Incisal

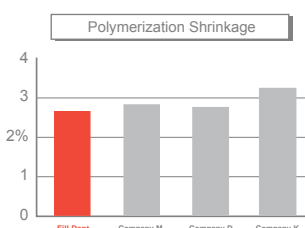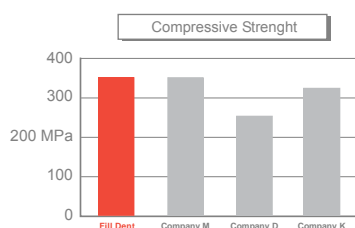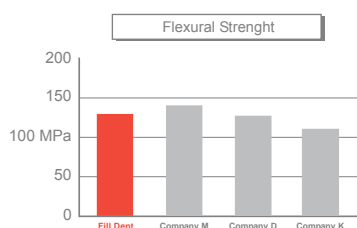

Syringes	0000100	Fill Dent Kit 5 Sir. + Bonding + Etch.
	0000101	Fill Dent Syr. Col. A1 4gr.
	0000102	Fill Dent Syr. Col. A2 4gr.
	0000103	Fill Dent Syr. Col. A3 4gr.
	0000104	Fill Dent Syr. Col. A3.5 4gr.
	0000105	Fill Dent Syr. Col. A4 4gr.
	0000106	Fill Dent Syr. Col. B1 4gr.
	0000107	Fill Dent Syr. Col. B2 4gr.
	0000108	Fill Dent Syr. Col. B3 4gr.
	0000109	Fill Dent Syr. Col.. C2 4gr.
	0000110	Fill Dent Syr. Col. C3 4gr.
	0000111	Fill Dent Syr. Col. D3 4gr.
	0000112	Fill Dent Syr. Col. A2 opaque 4gr.
	0000113	Fill Dent Syr. Col. A3 opaque 4gr.
	0000114	Fill Dent Syr. Col. EC 4gr.
	0000115	Fill Dent Syr. Col. T 4gr.
	0000116	Fill Dent Syr. Col. EI 4gr.

Capsule	0000117	Fill Dent Kit Cps
	0000118	Fill Dent Cps. Col. A1
	0000119	Fill Dent Cps. Col. A2
	0000120	Fill Dent Cps. Col. A3
	0000121	Fill Dent Cps. Col. A3,5
	0000122	Fill Dent Cps. Col. A4
	0000123	Fill Dent Cps. Col. B1
	0000124	Fill Dent Cps. Col. B2
	0000125	Fill Dent Cps. Col. B3
	0000126	Fill Dent Cps. Col. C2
	0000127	Fill Dent Cps. Col. C3
	0000128	Fill Dent Cps. Col. D3
	0000129	Fill Dent Cps. Col. A2 opaque
	0000130	Fill Dent Cps. Col. A3 opaque
	0000131	Fill Dent Cps. Col. EC
	0000132	Fill Dent Cps. Col. T
	0000133	Fill Dent Cps. Col. EI

# FILL DENT FLOW

LIGHT-CURED RADIOPAQUE FLOWABLE COMPOSITE RESIN



Fill Dent Flow is a light-cured radiopaque flowable composite resin, which has low shrinkage and good mechanical properties. Fill Dent Flow is designed to provide optimum flowability and handling for different types of restoration: direct restorations of Class V, Class III and small Class IV; as a liner for Class I and Class II restorations; repair of marginal or small defects in aesthetic indirect restoration; pit and fissure sealant; repair of temporary materials. It contains micro-hybrid fillers to provide excellent color adaptability, nature-like esthetics and high-gloss polishability.

## Indications

- Restorations of Class III, Class IV, Smaller Class V
- Base/liner in the Class I, Class II restorations
- Repair of resin, porcelain, and acrylic temporary materials
- Pit and fissure sealant
- Restoration of minimally invasive cavity preparations

## Packages

### Kit





- 2g x 4syringes (A2, A3, A3.5, B2)
- 40ea disposable tips (black)

### Refill

- 2g x 1syringe(A1, A2, A3, A3.5, A4, B2, B3, C2, A2O, A3O, I, A2L, A3L)
- 10ea disposable tips (black)



0000134	Fill Dent Flow Kit Syr 4x2gr. (A2-A3-A3.5-B2)
0000135	Fill Dent Flow Syr. Col. A1 2gr.
0000136	Fill Dent Flow Syr. Col. A2 2gr.
0000137	Fill Dent Flow Syr. Col. A3 2gr.
0000138	Fill Dent Flow Syr. Col. A3.5 2gr.
0000139	Fill Dent Flow Syr. Col. A4 2gr.
0000140	Fill Dent Flow Syr. Col. B2 2gr.
0000141	Fill Dent Flow Syr. Col. B3 2gr.
0000142	Fill Dent Flow Syr. Col. C2 2gr.
0000143	Fill Dent Flow Syr. Col. A2 opaque 2gr.
0000144	Fill Dent Flow Syr. Col. A3 opaque 2gr.
0000145	Fill Dent Flow Syr. Col. Incisal 2gr.
0000146	Fill Dent Flow Syr. Col. A2L 2gr. Low Viscosity
0000147	Fill Dent Flow Syr. Col. A3L 2gr. Low Viscosity

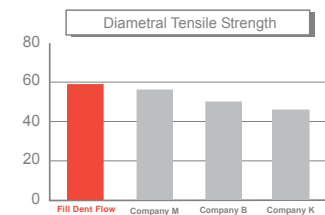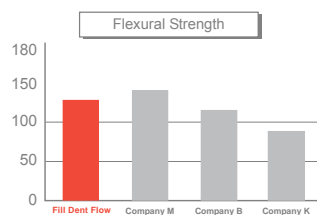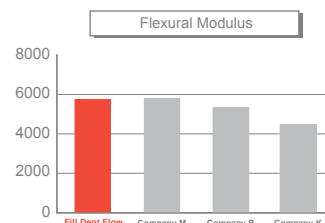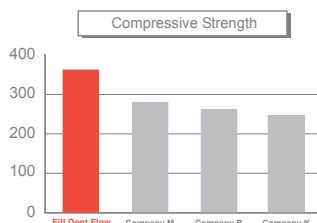

# ORTHO FLOW

LIGHT CURING ORTHODONTIC FLOWABLE COMPOSITE

Light curing composite for ceramic and metal brackets bonding, also indicated for retainers.

It is composed of Bis GMA/TEGDMA, with barium glass and silica. The content of inorganic filler (mean particle size is 0.01–2.5 mm) is 60% by weight.

## Advantages

- No additional mixing
- Can be easily and precisely applied on the bracket base
- Low viscosity
- Can be used for bonding orthodontic brackets without the use of a preliminary liquid resin on the etched enamel surface.
- Reduces working time

## Packages

Kit

4 syringes x 2g + 40 tips

Refill


1 syringe x 2g



00340144 | Ortho Flow Bond KIT 4 syringes x 2 gr.

00340174 | Ortho Flow Ref. Syr. 2gr.

*In vitro and in vivo studies confirm the advantages of the usage of Ortho Flow in orthodontic procedures.*

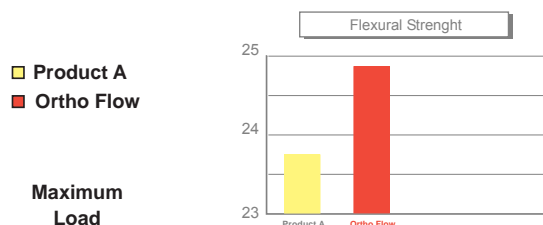

	Maximum Load
Product A	23,7125
Ortho Flow	24,9785

See Bibliography: Angle orthodontic vol. 75 n°3, 2005 - Angle orthodontic vol. 75 n°4, 2005



# ULTRA ONE BOND

ONE COMPONENT SELF-ETCHING LIGHT-CURED ADHESIVE



Ultra One Bond is 7th generation one-step bonding system with NANOFILLER technology.

Ultra One Bond is a one component self-etching light-cured adhesive designed to bond composite to dentin and enamel for direct and indirect restorations.

It uses NANOFILLER technology and shows a high adherence without separation. It also uses 4-META components with phosphate monomer.

## Packages:




Kit

- 5ml x 1 bottle
- Brush holder
- Brushes

Refill

- 5ml x 1 bottle

The film thickness of Ultra One Bond is less than 10 micrometer and it is homogeneous.

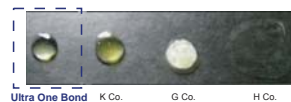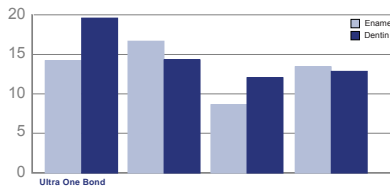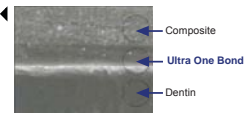

Using controlled surface coating technology of nano particles, there is no separation phenomenon during an evaporation of solvent.

0000162 | Ultra One Bond Kit

0000163 | Ultra One Bond Refill

# BASIC BOND

SINGLE COMPONENT BONDING AGENT



Basic Bond is the 5th generation single component bonding agent, light-activated, based on an ethanol formulation, designed to bond all classes of direct composite restorations to dentin, enamel, as well as procedures involving cast metals, composite, treated porcelain and set amalgam. Basic Bond allows indirect restoration procedures, due to its minimal and uniform film thickness.

## Indications:

- Bonding of direct composite
- Bonding to composite and set amalgam
- Bonding of indirect restorations-Porcelain
- Sealing hypersensitivity areas of teeth




## Packages

Refill

- 5ml x 1 bottle

## Microleakage

	Micro-leakage score	Standard Deviation
Basic Bond	1.1	0.9
S	1.0	0.8
P	1.6	1.3
O	1.0	1.2

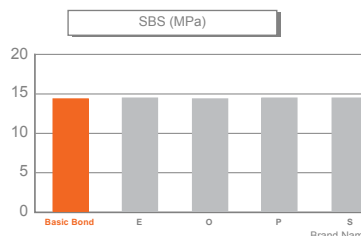

0000164 | Basic Bond Adhesive 5 ml



# FILL DENT SEAL

LIGHT CURED PIT & FISSURE SEALANT



0000157 Fill Dent Seal 4 syr.  
0000158 Fill Dent Seal syr. 1,2 ml

**Fill Dent Seal**, is a light-cured dental material designed for sealing the enamel pit and fissure of teeth for the prevention of caries. Thanks to specific nano-fillers, **Fill Dent Seal** has good properties like shear bond strength, micro-leakage and wear resistance. Its appropriate viscosity enables the sealant to penetrate the depths of pits. It doesn't need bonding

## Properties

Longer durability, low viscosity

Versatility, convenience, hygienic application in use  
Whiteness to be distinguishable from normal tooth structure.


## Packages

Kit

- Sealant 1.2ml x 2 syringes, 20 disposable tips - 2 syringes etchant, 20 disposable tips

Refill

Sealant 1.2ml x 1 syringe, 10 disposable tips

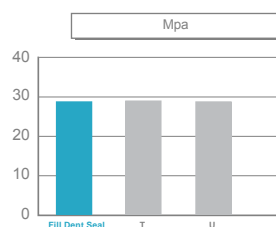

Materials	Micro-leakage score (standard deviation)
Fill Dent Seal	1.0000(1.1094)
T	0.9167(1.2795)
U	1.7333(1.4829)

# LIGHT TEMP

LIGHT-CURED TEMPORARY FILLING MATERIAL



**Light Temp** is a light-cured temporary filling material included nano-sized silver, indicated for temporary restorations. It has high elasticity and strength, it shows good marginal sealing properties. It can be easily used without any additional mixing and also can be easily removed without any debris in the cavity. Its color makes it distinguishable from normal tooth structure.




## Indications

- Temporary filling
- Inlay and onlay preparations

## Packages

- 3g x 3 syringes yellow (for anteriors)
- 3g x 3 syringes blue (for posteriors)

0000154 Light Temp YELLOW 3 syr. - 3gr.  
0000155 Light Temp BLUE 3 syr. - 3gr.

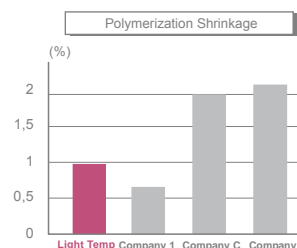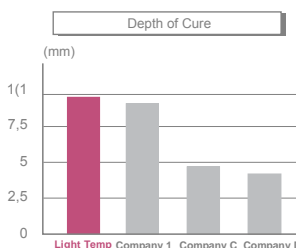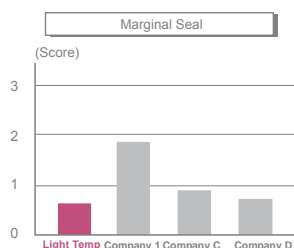

# FILL ETCH

ETCHING AGENT WITHOUT SILICA



**Fill Etch** is Phosphoric acid( $H_3PO_4$ ) without silica, with natural polymeric material.

A Semi-gel type etchant agent designed for easy etching work even in difficult points. Water-soluble and easy to remove; leaves no residues on the surface for a higher bonding strength on dentine and enamel.

## Characteristics:

- Biocompatible
- Leave no residues on the surface of etched tooth
- Easy and fast to remove
- 100% soluble

## Packages:

### Professional

- Etchant 5ml x 1 syringe + 10 disposable tips

### Professional plus


- Etchant 3ml x 2 syringes + 20 disposable tips

### Economic

- Etchant 3ml x 5 syringes + 20 disposable tips

### Bulk

- Etchant 30 ml x 1 syringe + 20 disposable tips

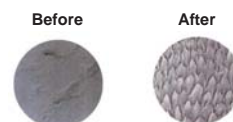

Microphotographs of premolar x 1,000

0000166	Fill Etch Professional 1 syr. x 5 ml
0000167	Fill Etch Professional Plus 2 syr. x 3 ml
0000168	Fill Etch Economic 5 syr. x 3 ml
0000169	Fill Etch Bulk 1 syr. x 30 ml

# OVER ETCH

NEW AND TOTAL ETCHING



**Over Etch** (35% Phosphoric acid) is a dental etching agent of viscous solution type, with thixotropic properties. **Over Etch** is designed to have a minimal cloggy effect during etching and to easily get all areas.



## Advantages

- Thixotropic property.
- Viscous solution type designed for easy etching work where you want.
- Extra fine tips suitable for canals, too.
- No residues on the surface of etched tooth

## Packages

- Etchant 1.2 ml x 4 syringes
- Etchant 30 ml x 1 syringe + dispenser



0000171	Over Etch Big Bulk 30ml x 1 syr.
0000172	Over Etch Kit 1,2ml x 4 syr.

# DUAL BUILD CORE

DUAL-CURED CORE BUILD-UP MATERIAL





**Dual Build Core** is a dual-cured (light+chemical) polymerization compounds. It is used as a replacement for dental caries and core build-up.

**Dual Build Core** is proper for the core build up due to its mechanical properties (compressive and flexural strength).

Thixotropic property of **Dual Build Core** is convenient for the dental treatment. It is easy to apply the treatment area where needed. And it is fully filled in a cavity with no empty space. NATURAL color available.

## Characteristics

- Compressive strength
- Thixotropic and long durable
- Auto mix type

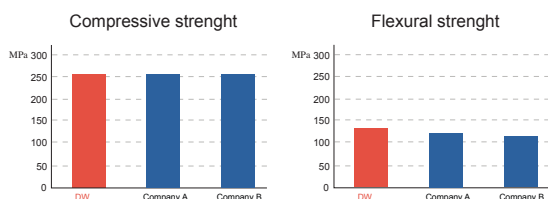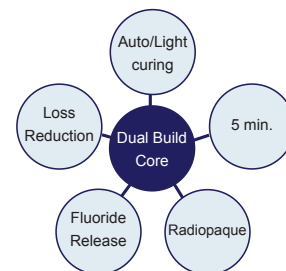

0000148	Dual Build Core N-Eco
000019	Dual Build Core N
0000150	Dual Build Core B-Eco
0000151	Dual Build Core B

# SELF BILD CORE

SELF-CURED CORE BUILD-UP MATERIAL



**Self Bild Core** is a self cure, fluoride releasing composite, used for core build-up.

Radiopaque, with excellent workability, smooth and creamy, easy to use thanks to the tips.

Short working time with low heat generation.

## Indications:


- Core Build-up
- Integral core and post procedure
- Post cementation
- Repairing and preventing root splitting

## Packages

- 50g x 1 Cartridge
- Automix tips

## Packages

- Base 10g + Catalyst 10g
- Mixing pad (50sheets), 15ea spatula
- 10ea disposable tip

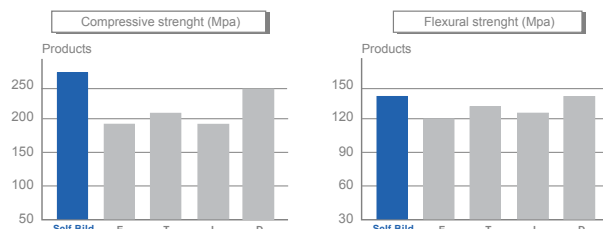

0000152	Self Bild Core 1 Cartridge (base 25g + 25g catalyst)
0000153	Self Bild Core (base 10g + 10g catalyst)

# ONE BUILD CEM

SELF-ADHESIVE UNIVERSAL RESIN CEMENT



0000159	One Build Cem W
0000160	One Build Cem U
0000161	One Build Cem O

**One Build Cem** is a dual-cured self-adhesive resin cement that requires no pre-treatment. It is easy to remove and has higher bond strength & radiopacity.

**One Build Cem** is available in Universal, White and Opaque.

## Packages:



- 9g x 1 automix syringe
- 10 intraoral tips
- 10 mixing tips

UNIFORM LUTING CEMENTATION

EASILY REMOVABLE

MULTI CEMENTATION

CONTINUOUS FLUORIDE ION RELEASE EFFECT

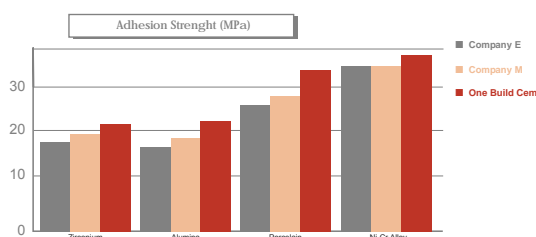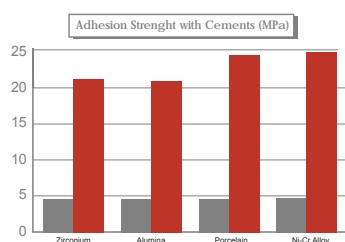

# CANAL PREP

CREAMY EDTA ROOT CANAL CLEANER



**Canal Prep** containing Ethylenediamine Tetraacetic acid (EDTA) as a chelating agent allows efficient cleaning and easy preparation of root canal. It contains urea-peroxide, whose bubbling effect efficiently lifts out the pulp debris in the canal.

## Indications:

- Shaping and cleansing of the root canal.
- Easy and fast removal of debris.
- Encouraging decalcification of dentin and antibacterial effect.
- Enables penetration of drugs into the root canal

## Package:

6 g. x 3 syringes

0000175 | Canal Prep E.D.T.A. 3 syr.x 6 gr.

# CANAL PASTE

CALCIUM HYDROXIDE PASTE



**Canal Paste** is water-based root canal filling material for infected root canal. It grants excellent flowability and accessibility. **Canal Paste** works on soft tissues and in the presence of liquids.

## Indications:

- Apexification
- Treatment for an infected root canal
- Temporary root canal filling
- Pulp capping of exposed pulp
- Disinfection after pulp extirpation treatment
- Pulp capping of vital pulpotomy

## Package:

- 2g x 1 syringe
- 20ea disposable tip / Extension holder
- 2g x 4 syringe

0000173	Canal Past Kit
0000174	Canal Past Refil



# FLUENT ELASTIC

IMPRESSION MATERIAL PUTTY



0000303 | Fluent Elastic- Violet 300+300 ml.

**Very high viscosity A-Silicone**

*for fixed prosthesis professional use in dental clinic.*

**Usage:** Recommended as special first impression material in case of undercuts for the two-stage technique in conjunction with a “wash” like Fluent Light or for the single-stage technique ( simultaneous) in conjunction with a “ wash” like Fluent Elastic

**Characteristic:**

- high elasticity and tear strenght
- high index of flowability
- doesn't compress the gingival mucosa
- excellent elastic recovery
- high resistance to deformation
- high dimensional stability
- no greasy surface
- violet colour

## Characteristics

Hardness	52 Shore
Iso 4823 - ADA 19	
Mixing Time	30 sec.
Time in mouth	150 sec.
Total Setting Time	3 min.
Mixing Ratio	1:1
Strain in compression	0,1 %
Elastic recovery	99,4 %

# FLUENT SOFT NORMAL

IMPRESSION MATERIAL PUTTY



0000302 | Fluent Soft Normal - Green 300+300 ml.

**Very high viscosity A-Silicone Bi-Component Super Fast Soft**  
*for fixed prosthesis professional use in dental clinic.*

**Usage:** Recommended as first impression material for all impression techniques : two-stage ( double impression) or single-stage (simultaneous) in conjunction with a “ wash” as Fluent Regular.

**Characteristic:**

- high index of flowability
- doesn't compress the gingival mucosa
- excellent elastic recovery
- high resistance to deformation
- high dimensional stability
- No greasy surface
- green colour

## Characteristics

Hardness	62 Shore
Iso 4823 - ADA 19	
Mixing Time	30 sec.
Time in mouth	150 sec.
Total Setting Time	3,5 min.
Mixing Ratio	1:1
Strain in compression	0,1 %
Elastic recovery	99,4 %

# FLUENT SOFT FAST

IMPRESSION MATERIAL PUTTY

Very high viscosity A-Silicone Bi-Component Super Fast Soft for fixed prosthesis professional use in dental clinic.

**Usage:** Recommended as first impression material in the two-stage technique ( double impression )

## Characteristic:

- rapid polymerization, shorter time in the oral cavity for a greater patient comfort
- high index of flowability
- doesn't compress the gingival mucosa
- excellent elastic recovery
- high resistance to deformation
- high dimensional stability
- no greasy surface
- yellow colour



0000301 | Fluent Soft Fast - Yellow 300+300 ml.

## Characteristics

Hardness	52 Shore
Iso 4823 - ADA 19	
Mixing Time	30 sec.
Time in mouth	90 sec.
Total Setting Time	2 min.
Mixing Ratio	1:1
Strain in compression	0,1 %
Elastic recovery	99,4 %

# FLUENT HARD

IMPRESSION MATERIAL PUTTY

Very high viscosity A-Silicone

for fixed prosthesis professional use in dental clinic.

**Usage:** Recommended as first impression material in the two-stage technique ( double impression )

## Characteristic:

- rapid polymerization, shorter time in the oral cavity for a greater patient comfort
- high index of flowability
- doesn't compress the gingival mucosa
- excellent elastic recovery
- high resistance to deformation
- high dimensional stability
- no greasy surface
- light blue colour



0000300 | Fluent Hard - Light Blue 300+300 ml.

## Characteristics

Hardness	62 Shore
Iso 4823 - ADA 19	
Mixing Time	30 sec.
Time in mouth	90 sec.
Total Setting Time	2 min.
Mixing Ratio	1:1
Strain in compression	0,1 %
Elastic recovery	99,4 %

# FLUENT LIGHT

IMPRESSION MATERIAL CARTRIDGE



0000304 | Fluent Light - Violet 2 cartridges 50gr.

Very light viscosity A-Silicone  
for professional use in dental clinic.

**Usage:** in the two-stage technique( double impression) in conjunction with a Putty of Extreme range

**Characteristic:**

- fast polymerization, shorter time in the oral cavity for a greater patient comfort
- high index of flowability
- excellent reproduction of details
- excellent elastic recovery
- high resistance to deformation
- high dimensional stability
- high tear strenght
- hydrophilic -thixotropic
- non-toxic -disinfectable
- violet colour

## Characteristics

Hardness	45 Shore
Iso 4823 - ADA 19	
Mixing Time	Automix
Time in mouth	150 sec.
Total Setting Time	2,5 min.
Strain in compression	0,1 %
Elastic recovery	99,8 %
Reproduction of details	2 micron

# FLUENT REGULAR

IMPRESSION MATERIAL CARTRIDGE



0000305 | Fluent Regular - Brown 2 cartridges 50gr.

Medium Viscosity hydrophilic addition silicone  
for professional use in dental clinic.

**Usage:** Recommended for precision impression as correction fluid suitable for the two-stage technique ( double impression) or as a monophase material in conjunction with a putty like Fluent Regular.

**Characteristic:**

- fast polymerization for best patient comfort
- high index of flowability
- excellent reproduction of the details
- excellent elastic recovery
- high resistance to deformation
- high dimensional stability
- high tear strenght
- hydrophilic -thixotropic
- non-toxic -disinfectable
- brown colour

## Characteristics

Hardness	56 Shore
Iso 4823 - ADA 19	
Mixing Time	Automix
Time in mouth	150 sec.
Total Setting Time	3 min.
Strain in compression	0,1 %
Elastic recovery	99,8 %
Reproduction of details	2 micron

# FLUENT OCCLUCAD

IMPRESSION MATERIAL CARTRIDGE

High viscosity addition silicone (polyvinylsiloxane) for a fast and accurate bite registration

**Usage:** Indicated for all intraoral registration, for registration of dental arches and for controls in centric relation. Compatible with CAD-CAM systems using optical, laser and contact reading.

## Characteristic:

- Fast polymerization for a major comfort for the patient
- High index of flowability
- Excellent reproduction of the details
- Fast < 60"
- High resistance to deformation
- High dimensional stability
- Hydrophilic -thixotropic
- Non-toxic -disinfectable
- Colour light blue



0000306 | Fluent Occlusal - Type I 2 cartridges 50gr.

## Characteristics

Hardness	95Shore
Iso 4823 - ADA 19	
Mixing Time	Automix
Time in mouth	20 sec.
Total Setting Time	60 sec.
Strain in compression	0,1 %
Elastic recovery	99,8 %
Reproduction of details	2 micron

# FLUENT CHECK

IMPRESSION MATERIAL CARTRIDGE

Low viscosity hydrophilic addition silicone for professional use in dental clinic.

**Usage:** Most suitable for disclosing areas of poor fit of prostheses. Recommended for checking pressure spots and the accuracy of fit of crowns and bridges.

## Characteristic:

- milky white opaque
- high index of flowability
- quick, easy accurate auto-mixing application
- minimum waste of material thanks to mini mixing tips
- in case of few elements, possibility of manual mixing
- optimized mixing time
- doesn't create thickness
- non-toxic
- white colour



0000307 | Fluent Check - White 2 cartridges 50gr.

## Characteristics

Hardness	80 Shore
Iso 4823 - ADA 19	
Mixing Time	Automix
Time in mouth	50 sec.
Total Setting Time	2 min.
Strain in compression	0,2 %
Elastic recovery	99,8 %
Reproduction of details	2 micron



# FLUENT CLEAR TIXO

Transparent silicon Automix **CLEAR**

for professional use in dental clinic and dental laboratory

IMPRESSION MATERIAL CARTRIDGE



**Usage:** Indicated in the realization of masks for the moulding of light-curing materials, suitable in free-hand technique in dentistry and dental laboratory, thanks to its high thixotropy.

## Characteristic:

- excellent thixotropy
- high index of fluency
- very high transparency level
- compact and smooth surfaces
- calibrated tear strenght
- high final hardness (70 Sh.A)
- dimensionally stable in time

[0000317](#) | Fluent Clear Tixo - White 2 cartridges 50gr.

# FLUENT IMPLA GUM

SEMI-RIGID TYPE



Addition vinylpolysiloxane for gingival reproduction.

**Usage:** Indicated for implant prosthetics for direct method realization of gingival reproduction in pink vinylpolysiloxane easy to trim, for the determination of the shape of the cervical contours and the emergence profile in fixed or mobile prosthesis

## Characteristic:

- trimming after 30 minutes from extrusion
- high index of flowability
- translucent colour which give it a natural appearance like real gingiva
- elasticity comparable to natural tissue
- balanced tear strenght
- high dimensional stability
- no greasy surface

## Characteristics

Hardness	70 Shore
Iso 4823 - ADA 19	
Mixing time	Automix
Working time	60 sec.
Time molding	10 min.
Time under drill	20 min.
Strain in compression	0,1 %
Elastic recovery	99,8 %
Reproduction of detail	2 micron

[0000319](#) | Fluent Impla 50 gr.

# DABLIU MIX

MIXING GUN



Mixing Guns for cartridge types. The gun fits all standard 50 ml cartridges.

[16/M3091](#) | Dabliu mix mixing gun 1:1 - 2:1  
[16/M3092](#) | Dabliu mix mixing gun 4:1 - 10:1

# DABLIU SIL

IMPRESSION MATERIAL



## Dabliu Sil

for professional use in dental clinic and dental laboratory.

**Usage:** indicated for technical masks, realization of gingival reproduction in pink vinylpolysiloxane. For mixed technique impression, to be used with **DABLIU SIL LIGHT** or **DABLIU SIL REGULAR**.

**0000308** | Dabliu Sil 62 Sh.A 1600 gr.

# DABLIU GUM

IMPRESSION MATERIAL



## Dabliu Gum

for professional use in dental clinic and dental laboratory.

**Usage:** indicated for technical masks, realization of gingival reproduction in pink vinylpolysiloxane. For mixed technique impression, to be used with **DABLIU SIL LIGHT** or **DABLIU SIL REGULAR**.

**0000309** | Dabliu Gum 52 Sh.A 1600 gr.

# DABLIU CATALYST

PASTE CATALYST FOR SILICONE RESINS

For technical use. Polydimethylsiloxane based catalyst for condensation silicone.



**0000312** | Dabliu Catalyst Paste 60 ml.

LIQUID CATALYST FOR SILICONE RESINS

Polydimethylsiloxane based catalyst for condensation silicone.

**Feature:** fluid ☐  
**Color :** clear



**0000313** | Dabliu Catalyst Liquid 30 ml.

# DABLIU SIL LIGHT

IMPRESSION MATERIAL



## Dabliu Sil Light

for professional use in dental clinic.

**Usage:** Polycondensation silicone. Indicated for dental impression as correction material. Precise reproduction of the finest details. No toxic release.

## Characteristics

Mixing Time	30 sec.
Time in mouth	180 sec.
Total Setting Time	210 sec.
Elastic recovery	99,4 %
Shore A	32/37

**0000310** | Dabliu Sil Light 150 gr.

# DABLIU SIL REGULAR

IMPRESSION MATERIAL



## Dabliu Sil Regular

for professional use in dental clinic.

**Usage:** Polycondensation silicone. Indicated for dental impression as correction material. Precise reproduction of the finest details. No toxic release.

## Characteristics

Mixing Time	30 sec.
Time in mouth	180 sec.
Total Setting Time	210 sec.
Elastic recovery	99,4 %
Shore A	40/45

**0000311** | Dabliu Sil Regular 150 gr.

# DABLIU ALGIN

HIGH PRECISION CHROMATIC ALGINATE



## Characteristics

Iso	1563
Type	Chromatic
Mixing time	35 sec.(pink)
Working time	90 sec.(white)
Setting time	120 sec.(blu)
Dosage	19 gr powder 40 ml water
Temperature	23 °C

### Notes:


Reported data are valid at 23°C. Temperature changes involve working and setting time changes as well.

Dabliu Algin is a high quality chromatic alginate for elastic high precision impressions. Indicated for models and fixed or removable prosthesis. Dabliu Algin reproduces details until 25 micron; lead- and cadmium-free. Tested with EN21563/ISO1563.

### Package:

- Bag x 450 g
- Box 25 bags

001025 | Chromatic Alginate 450 gr.

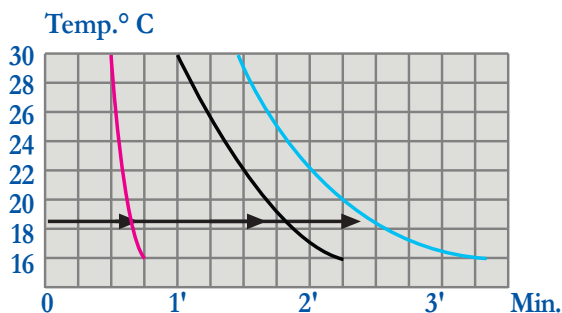

# DABLIU ALGIN GREEN

HIGH PRECISION NON-CHROMATIC ALGINATE



001032 | Non-Chromatic alginate 450 gr.

### Notes:

Reported data are valid at 23°C. Temperature changes involve working and setting time changes as well.

Dabliu Algin is a high quality alginate for elastic high precision impressions. Power-free, mint flavor, light green colour. Indicated for models and fixed or removable prosthesis. Lead- and cadmium-free; tested with EN21563/ISO1563.

### Package:

- Bag x 450 g
- Box 25 bags

## Characteristics

Iso	1563
Type	Non Chromatic(green)
Mixing time	35 sec.
Working time	80 sec.
Setting time	120 sec.
Dosage	19 gr powder 40 ml water
Temperature	23 °C

# ALGINBOY

ORTHODONTIC ALGINATE



High resistant orthodontic alginate for precise impression. Thanks to its stability, it is possible to braze the model even after 5 days ( if correctly stored).

### Characteristics

- Fruit flavour - Hypoallergenic - Rapid setting and pleasant taste
- Cadmium and lead free - Biocompatible - Chalk compatible
- Without powder for an hygienic use - Resistant to compression and breaking.

**\* Available in two different packages, bag or plastic drum**

001060 | Alginboy 1 drum 450gr.

001061 | Alginboy 1 bag 450gr.

# ENDO PERFECT

ENDODONTIC POINTS

Paper and Gutta Percha Points are hand-rolled sterilized endodontic points. Guaranteed consistency, size, shape and radiopacity in accordance with ISO6877, sterilization assurance level 10. Colored tip, packed in 28 mm boxes.

**Square Vial Box**



Note:

- Slide Box package available on request.
- Length mark available on request.

## - Absorbent Paper Point STANDARD SIZE (ISO/AEE)

Paper Points Endo Perfect 200pcs.		
	Size	Code
●	08	0000369
●	10	0000370
○	15	0000371
●	20	0000372
●	25	0000373
●	30	0000374
●	35	0000375
●	40	0000376
○	45	0000377
●	50	0000378
●	55	0000379
●	60	0000380
●	70	0000381
●	80	0000382
	Ass. 15/40	0000383
	Ass. 45/80	0000384

## - Absorbent paper point 5 pcs sterile blister (ISO)



Sterile Paper Points 180pcs.		
	Size	Code
●	08	0000477
●	10	0000478
○	15	0000479
●	20	0000480
●	25	0000481
●	30	0000482
●	35	0000483
●	40	0000484
○	45	0000485
●	50	0000486
●	55	0000487
●	60	0000488
●	70	0000489
●	80	0000490
	Ass. 15/40	0000491

## - Absorbent paper point NON STANDARD SIZE



Paper Points Endo Perfect 200pcs.		
	Size	Code
	XX-Fine	0000463
	X-Fine	0000464
	Fine	0000465
	Medium	0000466
	Coarse	0000467
	X-Coarse	0000468
	Assort	0000469





Note:

- Sterilize cell pack 180 pcs. available on request.
- Length mark available on request.

- Absorbent paper point TAPER SIZE; Taper 0.04 and Taper 0.06 (ISO)

Paper Points Endo Perfect 100pcs.		
Size	04%	06%
15	0000435	0000449
20	0000436	0000450
25	0000437	0000451
30	0000438	0000452
35	0000439	0000453
40	0000440	0000454
Ass. 15/40	0000441	0000455

### Square Vial Box



Note:

- Slide Box package available on request.
- Length mark available on request.

- Gutta Percha Points (ISO)

Gutta Percha Points 120pcs.	
Size	Code
08	0000222
10	0000223
15	0000224
20	0000225
25	0000226
30	0000227
35	0000228
40	0000229
45	0000230
50	0000231
55	0000232
60	0000233
70	0000234
80	0000235
Ass. 15/40	0000236
Ass. 45/80	0000237



- Gutta Percha Points NON STANDARD SIZE

Gutta Percha Points 100pcs.	
Size	Code
X-Fine	0000323
Fine-Fine	0000324
Medium-Fine	0000325
Fine	0000326
Fine-Medium	0000327
Medium	0000328
Medium-Large	0000330
Ass. XF - M	0000329

Gutta Percha Points 80pcs.	
Large	0000331
Gutta Percha Points 70pcs.	
X-Large	0000332



Note:

- Length mark available on request.


- Gutta Percha Points ISO TAPER SIZE; Taper 0.04 and Taper 0.06

Gutta Percha Points 60pcs.		
Size	04%	06%
15	0000288	0000302
20	0000289	0000303
25	0000290	0000304
30	0000291	0000305
35	0000292	0000306
40	0000293	0000307
Ass. 15/40	0000294	0000308



## DABLIU FILES NI-TI

NI-TI CANAL INSTRUMENTS

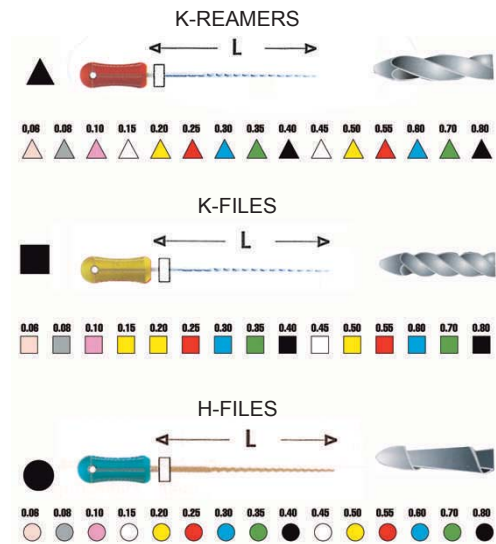

Ni-Ti instruments manufactured in an accurate process, in order to assure a total flexibility with shape memory. Their super elastic body guarantees the shape memory without altering qualities. ISO colour code.

**Package:**  
6 pcs. pack

21 mm	Dabliu Files NI-TI		
ISO	K-Reamers	K-File	H-File
08	144443800	144447800	144440800
10	144443801	144447801	144440801
15	144443802	144447802	144440802
20	144443803	144447803	144440803
25	144443804	144447804	144440804
30	144443805	144447805	144440805
35	144443806	144447806	144440806
40	144443807	144447807	144440807
45	144443808	144447808	144440808
50	144443809	144447809	144440809
55	144443810	144447810	144440810
60	144443811	144447811	144440811
70	144443812	144447812	144440812
80	144443813	144447813	144440813
Ass. 15/40	144443815	144447814	144440814
Ass. 45/80	144443814	144447815	144440815

25 mm	Dabliu Files NI-TI		
ISO	K-Reamers	K-File	H-File
08	144444800	144448781	144441791
10	144444801	144448782	144441792
15	144444802	144448792	144441802
20	144444803	144448802	144441812
25	144444804	144448812	144441822
30	144444805	144448822	144441832
35	144444806	144448832	144441842
40	144444807	144448842	144441852
45	144444808	144448843	144441853
50	144444809	144448844	144441854
55	144444810	144448845	144441855
60	144444811	144448846	144441856
70	144444812	144448847	144441857
80	144444813	144448848	144441858
Ass. 15/40	144444814	144450262	144450462
Ass. 45/80	144444815	144448849	144441859

28 mm	Dabliu Files NI-TI		
ISO	K-Reamers	K-File	H-File
08	144445800	144449800	144442800
10	144445801	144449801	144442801
15	144445802	144449802	144442802
20	144445803	144449803	144442803
25	144445804	144449804	144442804
30	144445805	144449805	144442805
35	144445806	144449806	144442806
40	144445807	144449807	144442807
45	144445808	144449808	144442808
50	144445809	144449809	144442809
55	144445810	144449810	144442810
60	144445811	144449811	144442811
70	144445812	144449812	144442812
80	144445813	144449813	144442813
Ass. 15/40	144445814	144449814	144442814
Ass. 45/80	144445815	144449815	144442815

31 mm	Dabliu Files NI-TI		
ISO	K-Reamers	K-File	H-File
08	144446800	144445100	144442900
10	144446801	144445101	144442901
15	144446802	144445102	144442902
20	144446803	144445103	144442903
25	144446804	144445104	144442904
30	144446805	144445105	144442905
35	144446806	144445106	144442906
40	144446807	144445107	144442907
45	144446808	144445108	144442908
50	144446809	144445109	144442909
55	144446810	144445110	144442910
60	144446811	144445111	144442911
70	144446812	144445112	144442912
80	144446813	144445113	144442913
Ass. 15/40	144446814	144445114	144442914
Ass. 45/80	144446815	144445115	144442915




## DABLIU FILES STEEL

STAINLESS STEEL CANAL INSTRUMENTS

Stainless steel instruments manufactured in a special stainless alloy steel which guarantees strength and total flexibility. Ergonomic handle for an excellent hold. Rubber stopper. ISO colour code.

Package: □

6 pcs. pack

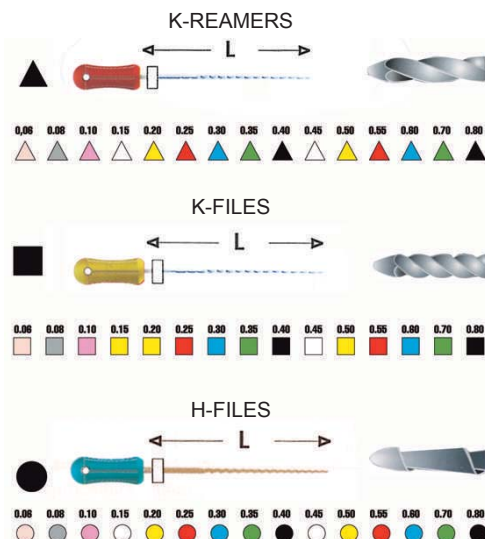

21 mm	Dabliu Files Steel		
ISO	K-Reamers	K-File	H-File
08	144517372	144518642	144515802
10	144517382	144518652	144515812
15	144517392	144518662	144516502
20	144517402	144518672	144516512
25	144517412	144518682	144516522
30	144517422	144518692	144516532
35	144517432	144518702	144516542
40	144517442	144518712	144516552
45	144517452	144518722	144516562
50	144517462	144518732	144516572
55	144517472	144518742	144516582
60	144517482	144518752	144516592
70	144516912	144515912	144516882
80	144516922	144515922	144516892
Ass. 15/40	144550022	144550222	144550422
Ass. 45/80	144550032	144550232	144550432

25 mm	Dabliu Files Steel		
ISO	K-Reamers	K-File	H-File
08	144517502	144518772	144515832
10	144517512	144518782	144515842
15	144517522	144518792	144511802
20	144517532	144518802	144511812
25	144517542	144518812	144511822
30	144517552	144518822	144511832
35	144517562	144518832	144511842
40	144517572	144518842	144511852
45	144517582	144518852	144511862
50	144517592	144518862	144511872
55	144517602	144518872	144511882
60	144517612	144518882	144511892
70	144516932	144515932	144511902
80	144516942	144515942	144511912
Ass. 15/40	144550062	144550262	144550462
Ass. 45/80	144550072	144550272	144550472

28 mm	Dabliu Files Steel		
ISO	K-Reamers	K-File	H-File
08	144517632	144518902	144515862
10	144517642	144519163	144515872
15	144519185	144519164	144519176
20	144519186	144519165	144519177
25	144519187	144519166	144519178
30	144519188	144519167	144519179
35	144519189	144519168	144519180
40	144519190	144519169	144519181
45	144519191	144519170	144519182
50	144519192	144519171	144519183
55	144519193	144519172	144519184
60	144519194	144519173	144512022
70	144519195	144519174	144512032
80	144519196	144519175	144512042
Ass. 15/40	144550102	144550302	144550502
Ass. 45/80	144550112	144550312	144550482

31 mm	Dabliu Files Steel		
ISO	K-Reamers	K-File	H-File
08	144517762	144519032	144515892
10	144517772	144519042	144515902
15	144517782	144519052	144515412
20	144517792	144519062	144515422
25	144517802	144519072	144515432
30	144517812	144519082	144515442
35	144517822	144519082	144515452
40	144517832	144519102	144515462
45	144517842	144519112	144515472
50	144517852	144519122	144515482
55	144517862	144519132	144515492
60	144517872	144519142	144515502
70	144517882	144519152	144515632
80	144517892	144519162	144515642
Ass. 15/40	144550142	144550342	144550542
Ass. 45/80	144550152	144550352	144550552

# PASTE CARRIER



Stainless Steel manufactured Lentulo dental paste carrier, ideal for injecting cements and medicaments into canals, with spring. ISO colour code. Autoclavable.

## Package:

- 4 pcs.

Length 21 - 25 - 29 mm

Ass. 25/40 - 25 - 30 - 35 - 40



Paste Carrier			
Size	21mm	25mm	29mm
Ass. 25/40	144550915	144550925	144550935
025	144513505	144512100	144514010
030	144513515	144512110	144514020
035	144513525	144512120	144514030
040	144513535	144512130	144514040


# BARBED BROACH



Stainless steel made, plastic handle with silicone, which grants a perfect hold. Antiseptic. Length: 21mm. ISO colour code Colours: violet, white, yellow, red, light blue, green, black.

## Packages:

15 blister packs x 10pcs.

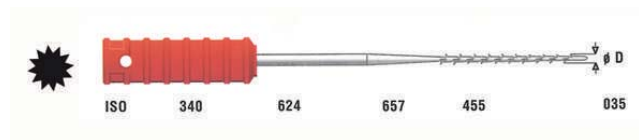

144515700S	Barbed Broach 10 pcs 15 grey
144516600S	Barbed Broach 10 pcs 20 violet
144516610S	Barbed Broach 10 pcs 25 white
144516620S	Barbed Broach 10 pcs 30 yellow
144516630S	Barbed Broach 10 pcs 35 red
144516640S	Barbed Broach 10 pcs 40 light blue
144516650S	Barbed Broach 10 pcs 50 green
144516660S	Barbed Broach 10 pcs 60 black


# ENDO ACCESSORIES



144510121 | Endo Clean



943540142 | Silicon Stops 100 pcs

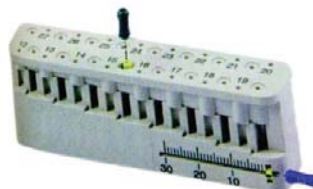

144510123 | Endo Bloc Gauge




144510122 | Endo Plastic Box



# DABLIU ENDO POST

ENDOCANAL POSTS

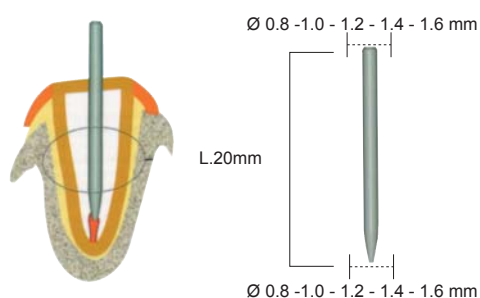

Dabliu Endo Post are carbon and glass fiber posts with high functional features. Realized with CNC modern technologies. Dabliu Endo Post, Conix Post and Crown Post are last generation endodontic posts and they can meet all operators' needs according to the treatment, anatomy and type of the root canal.

## Package:

10 pcs box.

## Sizes:

- Carbon: 0.8 - 1.0 - 1.2 - 1.4 - 1.6

- Glass: 0.8 - 1.0 - 1.2 - 1.4 - 1.6

## Kit

40 pcs. + 4 post drills

- Carbon: 0.8 - 1.0 - 1.2 - 1.4

- Glass : 0.8 - 1.0 - 1.2 - 1.4

0003801	Dabliu Endo Post CF 0,8 10 pcs
0003802	Dabliu Endo Post CF 1,0 10 pcs
0003803	Dabliu Endo Post CF 1,2 10 pcs
0003804	Dabliu Endo Post CF 1,4 10 pcs
0003805	Dabliu Endo Post CF 1,6 10 pcs
0003806	Dabliu Endo Post Kit CF 40 pcs 0.8-1.0-1.2-1.4 + 4 Drills


0003807	Dabliu Endo Post GF 0,8 10 pcs
0003808	Dabliu Endo Post GF 1,0 10 pcs
0003809	Dabliu Endo Post GF 1,2 10 pcs
0003810	Dabliu Endo Post GF 1,4 10 pcs
0003811	Dabliu Endo Post GF 1,6 10 pcs
0003812	Dabliu Endo Post Kit GF 40 pcs 0.8-1.0-1.2-1.4 + 4 Drills



Endo Post

# DABLIU CONIX POST

## ENDOCANAL POSTS

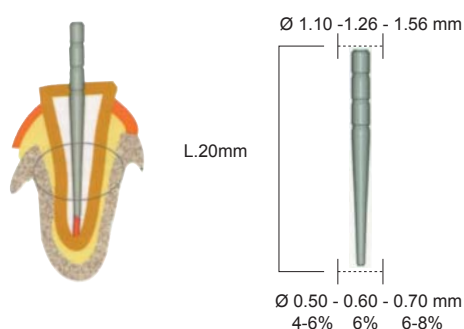

Dabliu Conix Posts are carbon or glass fiber posts, indicated for narrow roots and for small and medium restorations. Specifically realized for preparations with endodontic Ni-Ti instruments with retentive coronal part.

**Package:**  
10 pcs. box


**Sizes:**  
- Carbon: 0.50 - 0.60 - 0.70 mm  
- Glass : 0.50 - 0.60 - 0.70 - 0.80 mm



0003813	Dabliu Conix Post CF 0,50 10 pcs
0003814	Dabliu Conix Post CF 0,60 10 pcs
0003815	Dabliu Conix Post CF 0,70 10 pcs
0003816	Dabliu Conix Post GF 0,50 10 pcs
0003817	Dabliu Conix Post GF 0,60 10 pcs
0003818	Dabliu Conix Post GF 0,70 10 pcs
0003833	Dabliu Conix Post GF 0,80 10 pcs
0003836	Dabliu Conix Post GF 0,90 10 pcs

# DABLIU CROWN POST

## ENDOCANAL POSTS

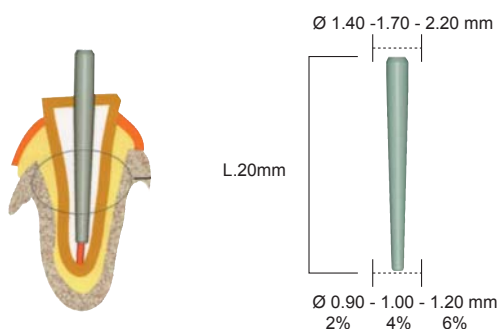

Dabliu Crown Posts is the ideal post for all kind of root canal. Indicated for preparations with the endodontic Crown Down technique.

**Package:**  
10 pcs box

**Sizes:**  
- Carbon: 0.90 - 1.00 - 1.20 mm  
- Glass : 0.90 - 1.00 - 1.20 mm



0003819	Dabliu Crown Post FC 0,90 2% 10 pcs
0003820	Dabliu Crown Post FC 1,00 4% 10 pcs
0003821	Dabliu Crown Post FC 1,20 6% 10 pcs
0003822	Dabliu Crown Post GF 0,90 2% 10 pcs
0003823	Dabliu Crown Post GF 1,00 4% 10 pcs
0003824	Dabliu Crown Post GF 1,20 6% 10 pcs

# SHORT POST R&R Micro

## ENDOCANAL POSTS



Short Post R&R Micro are high retention posts indicated for restorations, mainly to support and firmly fasten filling material to the tooth; they are also indicated for little core re-build.

### Package:


10 pcs. box

### Sizes:

Glass

dw	h	up
- 0.9	- 8.0	- 1.6 mm
- 0.9	- 9.0	- 1.6 mm
- 1.1	- 9.0	- 1.8 mm
- 0.9	- 1.0	- 1.6 mm
- 1.1	- 1.0	- 1.8 mm
- 1.1	- 1.1	- 1.8 mm

### Short Post R&R 890

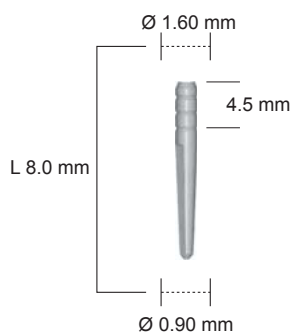

0003840 | Short Post R&R 890

0003841 | Short Post R&R 910

0003842 | Short Post R&R 911


0003843 | Short Post R&R 112

0003844 | Short Post R&R 113

0003845 | Short Post R&R 114


### Short Post R&R 910

### Short Post R&R 911



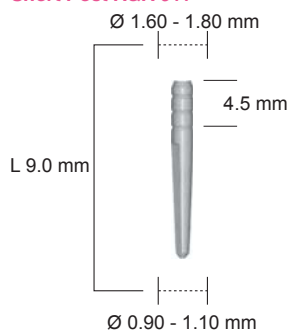


### Short Post R&R 112

### Short Post R&R 113

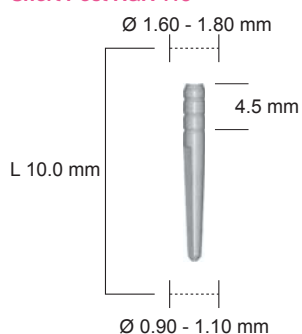

### Short Post R&R 114

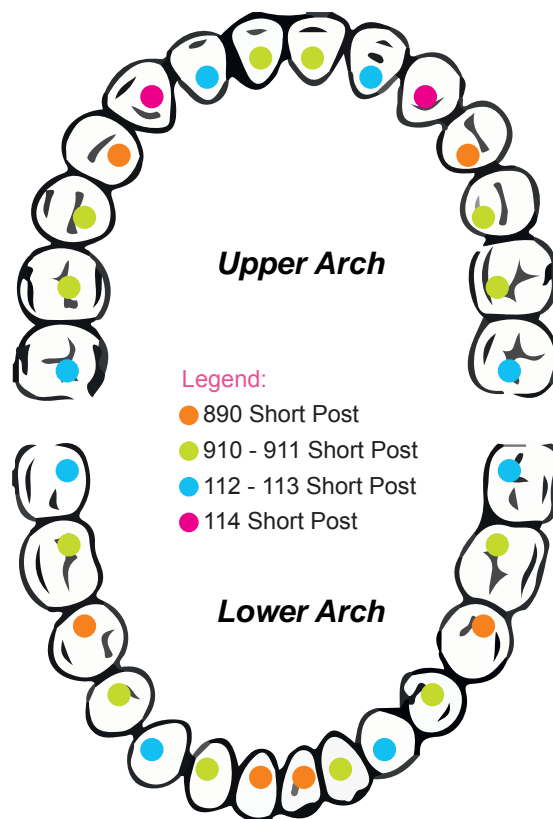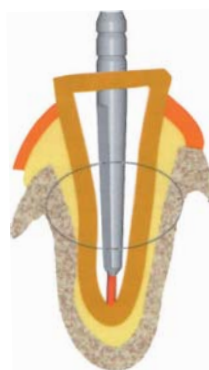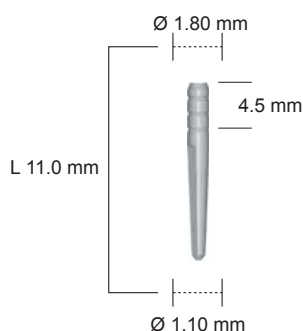

### Legend:

- 890 Short Post
- 910 - 911 Short Post
- 112 - 113 Short Post
- 114 Short Post

### Drills Guide

0003846	0003847
Short Post R&R 890	Short Post R&R 911
Short Post R&R 910	Short Post R&R 113
Short Post R&R 112	Short Post R&R 114

# DABLIU SMILE

BLEACHING GEL 38% (PEROXIDE HYDROGEN)

**Dabliu Smile 38%** is a peroxide hydrogen gel rich in potassium nitrate. Its bleaching action allows a 8 point chrome reduction. Its highly alkaline ph 9, got by mixing the two components, releases oxygen ensuring whitening action.

The potassium nitrate contained in Dabliu Smile, ensures low sensitivity and the best comfort for the patient. Its particular automix dispenser enables the separation between peroxide and activator, increasing self life.

## Package

- 2,5 ml x 3 syringes bleaching
- 3 ml x 1 syringe dam



Syringes

BEFORE



AFTER



0000015 | Dabliu Smile 38% Kit

# DABLIU DAM

LIGHT CURING LIQUID DAM

**Dabliu Dam** is a light curing liquid dam used to protect soft tissue and gums. Recommended for bleaching and restorative treatments. Blue color.

## Package:

- KIT
- 3.5 ml x 4 syringes + applicator tips
- BOX 30 KITS



Syringes



0000016 | Dabliu Dam Kit 4 syringes



# BRILLIANCE PEN

BLEACHING (CARBAMIDE PEROXIDE)

Brilliance is a carbamide peroxide based gel, mint flavored, used both for home and office bleaching. Very effective, well-tolerated, it causes no abrasion on the tooth surface. The presence of glycerin and of a desensitizer make the product less irritating than the hydrogen-peroxide based treatments. No sensitivity both during and after the bleaching treatment. Its mint flavor makes the product more agreeable and leaves a freshness sensation.

## HOW TO USE IT

- Proceed with an accurate teeth cleaning; use dental floss, if necessary. We suggest the use of Brilliance tooth paste, dental brush and dental floss;
- Carefully dry the tooth surface with cotton or tissue;
- Open the cap;
- Apply a thin layer on the anterior part of the teeth to be bleached and close the cap soon after;
- Keep the mouth open for about 30 seconds to let the gel dry;
- Do not eat or drink for about 30 minutes;
- After 30 minutes rinse the mouth with water and, if possible, wash your teeth.



0000018 | Brilliance Pen

## TREATMENT DURATION

Apply Brilliance Pen twice a day for one week, using the supplied brush. The product will be enough for the whole treatment. We suggest to daily wash your teeth with Brilliance whitening tooth paste that will help you to keep them white!



# BRILLIANCE

BLEACHING (CARBAMIDE PEROXIDE)

Brilliance is available in 4 different concentrations:

16% - 22% - 35% - 45%.

## 16% - 22% HOME BLEACHING

16% and 22% carbamide peroxide Brilliance are indicated for home bleaching through a mask, previously prepared by the dentist. The dentist gives the patient a syringe of Brilliance gel, too and patient carries on the treatment at home. The average period is 4-8 days, the mask must be kept in the mouth for 30-60 minutes.

## 35% OFFICE BLEACHING

35% carbamide peroxide Brilliance is indicated for office bleaching. It's very effective: it requires only one bleaching session!

## 45% OFFICE BLEACHING

45% for office bleaching can totally replace the hydrogen peroxide based treatments, as it's less aggressive and causes no abrasion on teeth surface!



0000002	Brilliance Bleaching Kit 16% 5ml x 4 syringes
0000003	Brilliance Bleaching Kit 22% 5ml x 4 syringes
0000014	Brilliance Bleaching Kit 35% 5ml x 3 syringes + DAM
0000008	Brilliance Bleaching Kit 45% 5ml x 3 syringes + DAM



# THERMOFORMING

MACHINE

To satisfy all professionals' needs, Horis dental has created two versions of this device: Silent and Aspy.

**SILENT** version represents the top model of Horis thermoforming machines range. It's automatic, silent and light, with an elegant design.

Silent is endowed with a last generation powerful and silent vacuum pump, which allows fast and precise realization of thermoformed devices.

The presence of a tray for metal pearls allows thermoforming process with gypsum models as well, without renouncing to optimal and professional results.

The plate for both square and round thermal discs, moves up and down thanks to an electric and silent motor.

All accessories included: basket, inox balls, electric cable and instruction manual.

**ASPY** version represents the economical model of Horis thermoforming line. It only differences from the Silent version for the fact that the plate for thermal discs moves manually.



00340080	Thermoforming Machine Aspy
00340700	Thermoforming Machine Silent

# TS-461

E-GENERATOR INTEGRATED LED HANDPIECE



**Description:** Self-illuminated, natural LED light. Can be used without fiber optic system.

Fits to universal dental units. Push button cartridge.

Autoclavable. One year guarantee (6 MONTHS FOR CARTRIDGE).

## TECHNICAL DATA:

OPERATION AIR PRESSURE: 0.25 MPa - 0.3 MPa

FREE LOAD SPEED: (280.000 TH - 350.000 SH) R.P.M.

SPRAY: Triple water spray

SUPPLY VOLTAGE: 3,0 - 3,8 V

GRIP POWER: 25 - 45 N

APPLICABLE BUR: (• 1.595- • 1.600 ) mm

NET WEIGHT: 88 g

SPEED OF GENERATOR: 18,000 - 22,000 R.P.M.

NOISE: not exceed 70 db

EXECUTION: according to ISO 7785-1:1999

LED LIGHT INTENSITY: • 1500 mcd

LED COLOR TEMPERATURE: 3500 K - 4000 K


LIFE OF LED BULB: • 5000 hours



1500056 | Turbo Speed Self Light Torque 1 spray LED 4 Holes

## TS-461 H

E-generator integrated LED handpiece

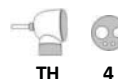

1500060 | Turbo Speed Self Light STD TS-461A LED 4 Holes

1500061 | Turbo Speed Self Light STD TS-461A LED 2 Holes

1500062 | Turbo Speed Self Light Torque TS-461A LED 4 Holes

1500063 | Turbo Speed Self Light Torque TS-461A LED 2 Holes

## TS-461 A

E-generator integrated LED handpiece



1500066 | Turbo Speed Self Light STD TS-461B LED 4 H. Quick Coupling

1500067 | Turbo Speed Self Light STD TS-461B LED 2 H. Quick Coupling

1500068 | Turbo Speed Self Light Torque TS-461B LED 4 H. Quick Coupling

1500069 | Turbo Speed Self Light Torque TS-461B LED 2 H. Quick Coupling

## TS-461 B

LED handpiece (integrated E-generator)  
Quick Connect Design



1500084 | Turbo Speed Self Light 45° TS-461C LED 4 Holes Surgical

1500085 | Turbo Speed Self Light 45° TS-461C LED 2 Holes Surgical (On Request)

1500057 | Turbo Speed 45° TS-461D 4 Holes Surgical

## TS-461 C

45 Degree surgical LED handpiece  
(E-generator integrated)






# TS-261

HANDPIECE WITH AND WITHOUT FIBER OPTIC

**Description:** Reliable performance via over 10,000 times press testing. Fiber optic illuminated handpiece with triple spray, brings comfortable feeling to patient. Use Hard fiber glass. Autoclavable. Easy to replace bulb. Unenclosed push button cartridge. Easy to change the bearing. Quick connector design, 360 degree swivel. One year guarantee (6 months for cartridge).

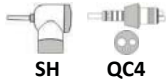


## TS-261

Classic without F.O.




### TECHNICAL DATA:

OPERATION AIR PRESSURE: 0.22 MPa - 0.3 MPa  
 FREE LOAD SPEED: (280.000 TH - 350.000 SH) R.P.M.  
 SPRAY: Triple Water Spray  
 SUPPLY VOLTAGE: 3,2 - 3,5 V  
 GRIP POWER: 25 - 45 N  
 APPLICABLE BUR: • 1.595 - • 1.600 mm  
 NET WEIGHT: 96 g  
 NOISE: not exceed 70 db  
 EXECUTION: according to ISO 7785-1:1999  
 LED LIGHT INTENSITY: • 200 Lux  
 LED COLOR TEMPERATURE: 3500 K-4000 K

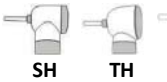

## TS-261 A-B

Classic with F.O.

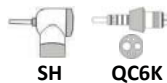

## TS-261 C

Classic with F.O.  
 Compatible with Kavo\*



1500070	Turbo Speed TS-261B Torque Quick Coup. F.O.
1500071	Turbo Speed TS-261A STD Quick Coup. F.O.
1500072	Turbo Speed TS-261 Torque Quick Coup. 4 H.
1500073	Turbo Speed TS-261C Quick Coup. Kavo* F.O.

# TS-221

HIGH SPEED AIR TURBINE DENTAL HANDPIECE FOR  
 INDIVIDUAL PATIENT

### Description:

1. Sterilized package, personal use, no cross infection
2. Maintenance free
3. Excellent design, fine appearance
4. 5 Hours working duration, be capable of treating 6-8 teeth
5. Quick connection (2 hole / 4 hole)



### TECHNICAL DATA:

OPERATION AIR PRESSURE: 0.22 MPa - 0.30 MPa  
 FREE LOAD SPEED: 350.000 R.P.M.  
 SPRAY: Single water spray  
 GRIP POWER: 25 - 45 N  
 APPLICABLE BUR: • 1.595 • 1.600 mm  
 NOISE: not exceed 70 dB  
 TORQUE: 0.08 N.CM  
 MATERIAL: handpiece handle-plastic  
 CARTRIDGE: copper & rubber  
 CHUCK TYPE: friction grip

1500091	Turbo Speed disposable handpiece TS-221 4 H. 10 pcs
1500093	Turbo Speed connector TS-221 4 Holes
1500094	Turbo Speed connector TS-221 2 Holes



## TS-421

RAINBOW HANDPIECE SERIES

Description: Standard &amp; torque head high speed handpiece

**GOLDEN**

4 Hole

**LIGHT BLUE**

4 Hole

**PINK**

4 Hole

**GREEN**

4 Hole

**BLUE**

4 Hole

**TECHNICAL DATA:**

OPERATION AIR PRESSURE: 0.20 MPa - 0.30 MPa

FREE LOAD SPEED: (350.000 - 400.000) R.P.M.

SPRAY: Single water spray

GRIP POWER: 25 - 45 N

APPLICABLE BUR: •1.595 ~ 1.600 mm

NET WEIGHT: 45 g


NOISE: not exceed 70 dB

EXECUTION: according to ISO 7785-1:1999

HANDLE MATERIAL: Aluminium alloys

CHUCK TYPE: push button

AVAILABLE COLOURS: golden, light blue, pink, green, blue.

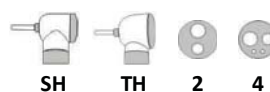

1500077	Turbo Speed	color TS-421 GOLD	4 Holes
1500078	Turbo Speed	color TS-421 PINK	4 Holes
1500079	Turbo Speed	color TS-421 LIGHT BLUE	4 Holes
1500080	Turbo Speed	color TS-421 BLUE	4 Holes
1500081	Turbo Speed	color TS-421 BLACK	4 Holes
1500082	Turbo Speed	color TS-421 GREEN	4 Holes

## LOW SPEED TS SERIES

HANDPIECE

**TS-9C**

Air motor (external water spray)



1500089 | Turbo Speed air motor TS-9C Ext. Spray 1:1

**TS-9D**

Air motor (inner water spray)



1500090 | Turbo Speed air motor TS-9D Int. Spray 1:1

**TECHNICAL DATA**

- TS 9C air motor

- 4 holes connector

- 1:1 low speed

- External water spray

SPEED: RPM 20.000

OPERATION AIR PRESSURE: 0,30~0,35 Mpa

NOISE: &lt; 70 dbm

**TECHNICAL DATA:**

OPERATION AIR PRESSURE: 0.30 MPa ~ 0.35 MPa

ROTATION SPEED: • 18.000 r/min

SPRAY: Inner water spray

NET WEIGHT: TS-7D: 85g

NOISE: • 70 dB


EXECUTION: according to ISO 7785-1:1999

# LOW SPEED TS SERIES

## HANDPIECE

### TS-75 / TS-77

E - generator integrated low speed handpiece (inner water spray)

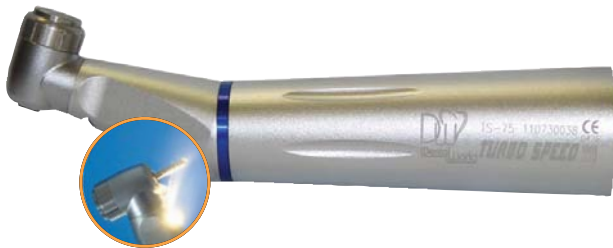

**1500086** | Turbo Speed Contra angle TS-75 Int. Spray 1:1 LED

**1500053** | Turbo Speed Contra angle TS-77 Int. Spray 1:1

#### TECHNICAL DATA:

OPERATION AIR PRESSURE: 0.30 MPa ~ 0.35 MPa

ROTATION SPEED: • 18.000 r/min

SPRAY: Inner water spray

SUPPLY VOLTAGE: 3.0 - 3.8 V

APPLICABLE BUR: TS-75 • 2.345 ~ 2.355mm

TS-76 • 1.595 ~ 1.600mm

NET WEIGHT: TS-75 / TS-77: 50g

NOISE: • 70 dB

EXECUTION: according to ISO 7785-1:1999

LED LIGHT INTENSITY: • 1500 mcd

LED COLOR TEMPERATURE: 3500 K - 4000 K

LIFE OF LED BULB: • 5000 hours

CHUCK TYPE: push button

### TS-73

Push button contra angle (high speed dental bur)



**1500087** | Turbo Speed Contra angle TS-73 Ext. Spray 1:1 CA

#### TECHNICAL DATA:

OPERATION AIR PRESSURE: 0.30 MPa ~ 0.35 MPa

ROTATION SPEED: • 18.000 r/min

SPRAY: External water spray

NET WEIGHT: TS-73 / TS-74: 42g

NOISE: • 70 dB

EXECUTION: according to ISO 7785-1:1999

### TS-79

Low Speed reduction



**1500058** | Turbo Speed Contra angle TS-79 CA 20:1

#### TECHNICAL DATA:

OPERATION AIR PRESSURE: 0.30 MPa ~ 0.35 MPa

CHUCK TYPE: push button

SPRAY: External irrigation system

NET WEIGHT: TS-79: 40g

NOISE: • 70 dB

AVAILABLE: according with ISO standard

COMPATIBLE: with all E-type motors

### TS-8B / TS-8C

Straight handpiece



**1500088** | Turbo Speed Straight Handpiece TS-8B Ext. Spray 1:1

#### TECHNICAL DATA:

Low speed handpiece 1:1 Ratio

OPERATION AIR PRESSURE: 0.30 MPa ~ 0.35 MPa

ROTATION SPEED: • 18.000 r/min

SPRAY: External water spray

NET WEIGHT: TS-8B: 44g

NOISE: • 70 dB

EXECUTION: according to ISO 7785-1:1999



**1500052** | Turbo Speed Straight Handpiece TS-8C Int. Spray 1:1

#### TECHNICAL DATA:

Low speed handpiece 1:1 Ratio

OPERATION AIR PRESSURE: 0.30 MPa ~ 0.35 MPa

ROTATION SPEED: • 18.000 r/min

SPRAY: Internal water spray

NET WEIGHT: TS-8C: 44g

NOISE: • 70 dB

EXECUTION: according to ISO 7785-1:1999

\* All Contra angle available for FG Burs

# SPORIDOX

DISINFECTANT



**Sporidox** is a decontaminant and disinfectant detergent solution based on active oxygen and peracetic acid, with a sporicidal, tuberculocidal, virucidal, fungicidal action. Used for decontamination, cleaning and disinfection of medical and surgical devices.

**Packages:**

- Jar 1000g
- BOX 6 jars

**Use:**

Decontamination and cleaning of medical and surgical devices (D.L. Sept 28, 1990). High level disinfection with contemporary cleaning of medical and surgical devices.

00340023 | Sporidox IMS

# SPOREXIN plus DETERSONIC

DISINFECTION BY ULTRASOUND



**Sporexin Plus Detersonic** is a disinfectant and detergent solution based on polyphenols, efficient as tuberculocidal, bactericidal, virucidal, fungicidal, included Mycobacterium tuberculosis, broad-spectrum antimicrobial activity of polyphenols is associated to detergent capacity. It keeps its action in presence of organic substances. Used to clean and disinfect medical devices, before sterilization. Ultrasounds use reduces cleaning and disinfection time.

**Packages:**

- Bottle 1000 ml
- BOX 12 bottles

**Use:**

Cleaning and disinfection of medical devices, before sterilization phase. Immersion of sharp and pointed instruments after their use, according to D.L. Sept. 28, 1990. Ultrasounds use reduces cleaning and disinfection time.

54/0721507 | Sporexin plus detersonic 1000ml IMS

# PROTEOZIM plus 400

DISINFECTION BY ULTRASOUND



Fast action enzymatic solution against blood, proteins and fats. Diluted 1:400 (2,5 ml in 1 lt water), it's active in 10 minutes in open trays; in 5 minutes in ultrasounds basin and automatic washer disinfectors.

**Packages:**

- Bottle 1000 ml
- Box 12 bottles


**Use:**

It's indicated for washing and cleaning of instruments, surgical and medical devices. It can be used for manual washing of endoscopes and glassware, in ultrasounds basin.

54/0720028 | Proteozim Plus 400 1lt

# PROTEAZONE OD

PRE-STERILIZATION

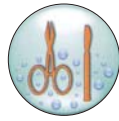

**Proteazone OD** is a concentrated decontaminant solution with balanced formula, synergized with Adazone for dentistry devices. Myco/bactericidal, virucidal (HIV, HBV, HCV) and fungicidal specific for microbial biofilm. Active in 5 minutes. Indicated for washing, disinfection of instruments, equipments, medical devices. It can be used in open trays and in ultrasound basin.

## Packages:

- 1000 ml
- Box 12 pcs.

## Use:

Washing, disinfection of instruments, equipments, medical devices. It can be used in open trays and in ultrasound basin. Dilute 1% (10ml in 1 l water) for activity in 10 minutes. Dilute 2% (20 ml in 1 l water) for activity in 5 minutes.

54/0720027 | Proteazone OD 1LT.

# ADASPOR READY

STERILIZATION



**Adaspor Ready** is a cold sterilizing solution, to be activated, based on Adazone and peracetic acid (PAA).

In 5 minutes it guarantees a sporicidal and high level disinfection. Sterilization in 10 minutes. Once activated, the solution can be used for 12 days. Ready-to-use solution for cold sterilization and / or high level disinfection medical devices. Particularly recommended in endoscopy centres, because of its fast action.

## Packages:

- Box with solution and activator (Sol. A+ Sol. B) 1000 ml.

## Use:

Cold sterilization and/or high level disinfection of medical devices.

500140 | Adaspor Ready

# SPOREXIN plus OD

STERILIZING AND HIGH LEVEL DISINFECTANT SOLUTION



**Sporexin Plus OD** is a concentrated solution for sterilization and high level disinfection. Glutaraldehyde and polyphenols based solution, it exercises a sporicidal, tuberculocidal, bactericidal, virucidal, fungicidal action, included Mycobacterium tuberculosis. Sterilization in 10 hours; high level disinfection in 10 minutes. It can be used in open trays.

## Packages:

- 4 bottles (250 ml)

## Use:

Cold sterilization and high level disinfection of dentistry devices. It can be used in open trays and in ultrasound basin. 1:10 dilution sterilizes in 10 hours, disinfects at high level in 10 minutes. Activate solution, kept in cover bin, lasts 30 days. Low concentration of glutaraldehyde.

54/0720026 | Sporexin plus OD IMS CF. 4x250ml.



## SPOREXIN<sup>plus</sup> DS

SPRAY FOR SURFACES



Sporexin Plus DS is a polyphenols, water based spray glutaraldehyde, with a spray dispenser. It has a tuberculocidal, bactericidal (included Mycobacterium tuberculosis) action. Sporexin plus is active also in presence of large amount of organic substances the solution presents a long residual bacteriostatic activity, if not removed from surfaces on which it is applied. Ready to use solution for intermediate level disinfection of operating theatre devices and dental devices and furniture.

### Packages:

- Case with 5 bottles, 500 ml and 2 spray closures
- Box 4 cases

### Use:

Work surfaces and benches, operating theatre devices such as lights, trolleys, operating tables, dental devices, furniture and instrument trays of dental suites. Ready to use.

[54/0721506](#) | Sporexin plus DS IMS pz 5x500ml.

## SPOREXIN<sup>plus</sup> DS towels

TOWELS FOR SURFACES



Sporexin Plus DS towels is a polyphenols solution without glutaraldehyde. It has a tuberculocidal, bactericidal, virucidal, fungicidal (included Mycobacterium tuberculosis) action. It is active also in presence of large amount of organic substances on which it is applied. Ready-to-use towels for intermediate level disinfection of work surfaces and benches: operating theatre devices such as lights, trolleys, operating tables hospital aisle beds furnitures.

### Packages:

- Dispenser with 120 towels
- Box 6 dispensers

### Use:

Work surfaces and benches, operating theatre devices such as lights, trolleys, operating tables, dental devices, furniture and instrument trays of dental suites.

[54/0721509](#) | Sporexin Plus DS towels

## ISACLEAN

SPRAY FOR SURFACES



Based on the latest generation patented molecule called ISAZONE. Rapid and efficient cleaning and disinfectant solution, ready to use. No carcinogen and toxic substances, easy to use, it is active in 10 minutes. The presence of Isazone, combined with benzalkonium chloride and chlorhexidine digluconate, increases its efficiency (after 2 minutes of activity).


### Packages:

- Case with 500ml x 5 bottles, 2 spray closures

[550020](#) | Isaclean Spray 5 x 500ml

## SPOREXIN<sup>plus</sup> VACUUM

DISINFECTANT FOR SOLUTION DEVICES

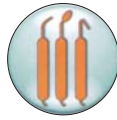

Sporexin Plus Vacuum is a polyphenols solution, specific for the disinfection of aspiration circuits with antifoam action. It reduces chalky depot formation with its tuberculocidal, bactericidal, virucidal, fungicidal (included the *Mycobacterium tuberculosis*) action. Recommended use: dissolved 6% for disinfection of suction circuits and dental surgery.

**Packages:**

- Tank 5000 ml

**Use:**

Indicated for dental suction devices disinfection.  
Dilute 6% (60ml in 1 l water)

54/0721504 | Sporexin vacuum 5LT.

## GERMICIDIN

DISINFECTANT, DETERGENT, CONCENTRATED SOLUTION

N-alkyl-dimethyl-benzyl-ammonium chloride 10%, EDTA 0.25% and inert ingredients. Indicated for floors and surfaces in general.



**Packages:**

Tank 5000 ml

**Use:**

Dilute the solution 1:250 (10ml in 2,5 l water)

54/0720024 | Germicidin Tank 5lt

## BLUSTERIL ALCOHOLIC

DISINFECTANT, CONSERVANT

Blusteril Alcoholic is an alcoholic solution based on benzalkonium chloride with antioxidants. It is a ready-made disinfectant. High level bactericidal, virucidal, fungicidal solution for medical devices. It is indicated for the disinfection and aseptic conservation of medical devices.



**Packages:**

- Bottle 1000 ml
- BOX 12 bottles

**Use:**

Disinfection and aseptic conservation of medical devices. Keep the solution in covered containers.

505025 | Blusteril alcoholic 1LT

# BENODENT

MOUTH WASH



**Benodent CLX 0.2** is indicated before and after extraction, small surgeries, extractions, antiplaque intensive and maintenance treatment. Benodent intensive use has a continuative antiplaque action; chlorhexidine keeps a long term protective action on mucosa, gingiva and enamel.

**Benodent FLU 0.02** is a mouthwash with 0.02 chlorhexidine and 0.07 sodium fluoride, specifically indicated for caries and gingivitis prevention, aids to reduce bacterial plaque formation.

## Packages:

- Bottle 1000 ml  
BOX 8 pcs.
- Bottle 500 ml  
BOX 6 pcs
- Single-dose 15 ml

3012	Benodent CLX 0.20% bottle 500 ml
0821508	Benodent CLX 0.20% 2 bottles 1000ml
3015	Benodent CLX 0.20% 10 pcs / pack 15 ml
3022	Benodent CLX 0.02% and FLU bottle 500 ml
0821509	Benodent CLX 0.02% and FLU 2 bottles 1000 ml
3025	Benodent CLX 0.02% and FLU 10 pcs / pack 15 ml

# SPOREXIN

HAND SOAP



**Sporexin hands soap** is an antiseptic detergent for hands, based on poly phenols. It has a tuberculocidal, bactericidal, virucidal and fungicidal (included Mycobacterium tuberculosis) solution. It is indicated for the presurgical antisepsis of hands and forearm.

## Packages:

- Bottle 500 ml  
BOX 24 pcs.
- Tank 5000 ml  
BOX 4 tanks

54/0721508	Sporexin hands 500 ml Bottle
54/0721505	Sporexin hands 5000 ml Tank

# DISINFECTION TRAY

Rectangular trays indicated for sterilizing and disinfecting dental devices, with internal basket and covering top.



00340900	Disinfection Tray 1 lt.
00340041	Disinfection Tray 3 lt.

# HAND SANITIZER GEL

Horisan Hand Sanitizer Gel is an antiseptic gel specifically designed to thoroughly sanitize the skin of the hands, removing 99,9% micro-organisms in every circumstance, without water.

## Use:

Horisan Hand Sanitizer Gel ensures protection thanks to the alcohol contained that, with its special formulation, protects and cares skin. It is practical in any situation, where you want, to clean or sanitize hands on public transport, in hospital, in public toilets, after handling money, phones, etc. It's particularly indicated for doctors and medical assistants.

Pour a small amount directly on your hands and rub until completely dry. Horisan Hand Sanitizer Gel is gently scented and doesn't leave your hands sticky.



0001005 | Horisan hand gel 500ml



01010 | Horisan hand gel 80ml



0001001 | Horisan hand gel 1000ml

## Cautions:

External use- Inflammable material

Avoid the contact with eyes. In case contact occurs, rinse with abundant water.

Keep the container closed, far from heat and from children. Do not drink.

Not to disperse the empty container in the atmosphere.

Ingredients: Alcohol denat, Aqua, Diisopropyl adipate, Glycerin, Parfum, Propylene cicol, Carbomer, Aminomethyl propanol

## Packages:

80 ml -500ml -1000ml.





Ortho Bite - Ortho Box



## ORTHO BOX & ORTHO BITE

Ortho box and Ortho bite are original, functional, easy containers to carry out orthodontic plaques, bite or functional appliances. Manufactured with no toxic material, they guarantee hardness and durability. Available in a wide range of beautiful colours. Two sizes: 40 mm and 30 mm.

### Ortho Bite

Blue	00340235
White	00340230
Violet	00340234
Orange	00340227
Fuxia Neon	00340228
Yellow Neon	00340233
Fuxia	00340232
Green	00340231
Red	00340236
Light Blue	00340229
Yellow	00340237
Black	00340239
<b>1 piece</b>	



### Ortho Box

Violet	00340210
Red	00340205
Green	00340202
White	00340201
Dark Yellow	00340209
Light Blue	00340207
Yellow Neon	00340203
Fuxia Neon	00340204
Light Yellow	00340215
Blue	00340206
Fuxia	00340211
Silver	00340208
Orange	00340199
<b>1 piece</b>	

## PROSTHESES CASE

Orthodontic basin and prosthesis case ideal for cleaning all type of prosthesis and orthodontic appliance. It can contain superior prostheses and inferior one. Filter for impurity endowed. Available in a wide range of beautiful colors.



### Prostheses Case

White	00771469
Orange	00771470
Blue	00771472
Sky Blue	00771472A
Light Blue	00771472C
Fuxia	00771472F
Yellow	00771472G
Red	00771472R
Violet	00771472V
Green	00771472VR
Fuxia Neon	00771472FN
Carton packed	00771472BOX
<b>1 piece</b>	

# DABLIU BRUSH

Micro disposable applicators made of not absorbent fiber, bendable, 10 cm long, available in round tubes or in square version, suitable for dispenser.

Size: Regular, Fine, Super Fine.

Package:

100pcs x 4 tubes.



00100101	Dabliu brush 4x100 Regular
00100102	Dabliu brush 4x100 Fine
00100103	Dabliu brush 4x100 Super Fine
00100105	Dabliu brush Dispenser 4x100 Regular
00100106	Dabliu brush Dispenser 4x100 Fine
00100107	Dabliu brush Dispenser 4x100 Super Fine
00100108	Dabliu brush Dispenser

# DABLIU ULTRA BRUSH

DISPOSABLE BRUSH APPLICATORS

Characteristics:

- Disposable
- Bendable
- Microfiber

Package:

Dabliu Ultra Brush, 100 pcs.



00100109	Dabliu ultra brush 100 pcs
----------	----------------------------

# DABLIU ENDO BRUSH

ENDODONTIC MICRO APPLICATORS

Characteristics:

- Disposable
- Bendable
- Microfiber

Package:

Dabliu Endo Brush, 100 pcs.



00100110	Dabliu endo brush regular 100 pcs
00100112	Dabliu endo brush fine 100 pcs

# MICRO-APPLICATORS

MICRO APPLICATORS FOR BRUSH HANDLE



## Characteristics:

- Disposable
- Brush handle autoclavable

## Package:

Brush handle 1 pc

## Package:

Micro applicators tips 100 pcs

00100111 | Micro - applicators 100 pcs

0034117 | Brush handle 1 pc

# VISOR GALAXY

Visor shield with removable screen realized with plastic support easily adaptable to heads with three screens supplied. The screens have a perfect visibility, then do not tarnish and they are cleanable with a soft cloth.

## Package:

BOX 12 pcs.



## Color shade



3P1011	Visor Galaxy White
3P1011F	Visor Galaxy Fuxia
3P1011G	Visor Galaxy Yellow
3P1011V	Visor Galaxy Violet
3P1011FN	Visor Galaxy Fuxia Neon
3P1012	Visor Galaxy Green
3P1013	Visor Galaxy Blue
3P10112A	Visor Galaxy Orange
3P1021	Screen Refill Visor Galaxy 3 pcs

# PROTECTIVE GLASSES

## Indications:

- Clear, to protect from liquid and squirts
- Orange, to protect from blue light radiations.

## Characteristics:

- Anti-scratch
- Anti-fog

## Package:

Protective glasses 1pc.



00400012 | Protective Glasses Orange

00400018 | Protective Glasses Clear

# PROTECT PLUS

DENTAL BIBS WITH POCKET



## Dispenser



Protect Plus thanks to their 2-ply cellulose and polyethylene, grant excellent resistance and absorbance.

Easy to wear and to use, they are individually pre-cut and endowed of a pocket for patient's comfort.

Packed in a practical dispenser to avoid contamination and contact with external agents.

Size: □

68 cm x 48 cm

Package:

100 pcs dispenser

### Plus resistance

2-ply cellulose and polyethylene



### Plus absorbent

high quality cellulose and pocket



### Plus comfort

size 68x48cm and easy to use



### Plus waterproof

thanks to polyethylene layer



00870001 | Protect Plus 68 cm x 48 cm 100 pcs

# DABLIU TIPS

APPLICATOR NEEDLES



00340158



00340112



00340111

Needle tips for syringes with Luer-Lock connection. □

1- BLUE Ga. 25 Ø 0,5 mm for etching material

2- Black Ga. 20 Ø 0.92 mm for flow material

3- Orange with piston Ga. 19 Ø 0.88 mm for cements

4- Needle for EDTA Ga. 20 Ø 0.98 mm

5-6 Ultrafine plastic needles for cements



00340137



00340119M



00340119

00340111	Dabliu Tips Blue Refill 20 pcs
00340111A	Dabliu Tips Blue Refill 100 pcs
00340112	Dabliu Tips Black Refill 20 pcs
00340112A	Dabliu Tips Black Refill 100 pcs
00340158	Dabliu Tips Violet Refill 20 pcs
00340137	Dabliu Tips Orange w/piston Refill 100 pcs
00340119	Dabliu Tips Plastic 100 pcs
00340119M	Dabliu Tips Plastic lenght mark 100 pcs



# DISPENSER FOR COTTON

Handy and hygienic dispenser for cotton rolls. Its new mechanism permits through a light pressure, the single dispensation of rolls.




00110012	Dispenser for cotton WHITE
00110017	Dispenser for cotton BLUE
00110018	Dispenser for cotton GREEN
00110019	Dispenser for cotton YELLOW

# CRISTAL TOOTH


Dental jewel in crystal, available in 9 gradation of colour, diameter 1.8 mm. Mechanical retention on base to obtain a better adhesion. Colour available:

Gold - Iris - Light Blue - Red - Pink - Black - Green - Violet - Crystal - Assorted

Package:  
10 pcs

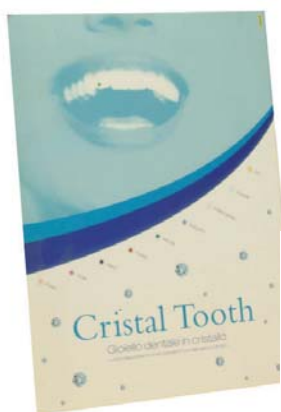

Display rack

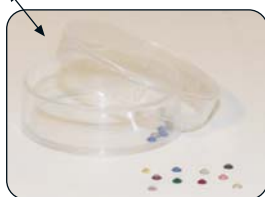

- Gold    ● Iris    ● Light Blue    ● Red    ● Pink  
● Black    ● Green    ● Violet    ● Crystal    ● Assorted

07001	Crystal Tooth 10 pcs ASSORTED
07004	Crystal Tooth 10 pcs CRYSTAL
07005	Crystal Tooth 10 pcs IRIS
07006	Crystal Tooth 10 pcs GOLD
07007	Crystal Tooth 10 pcs LIGHT BLUE
07008	Crystal Tooth 10 pcs GREEN
07009	Crystal Tooth 10 pcs RED
07010	Crystal Tooth 10 pcs BLACK
07011	Crystal Tooth 10 pcs VIOLET
07012	Crystal Tooth 10 pcs PINK

# APPLYTRAY

FLUORIDE GEL TRAYS

Soft fluoride gel tray for professional dental purpose.

Package:

200 pieces (100 upper + 100 lower)

- Small (white)
- Medium (yellow)
- Large (blue)



00340146	Eazy Tray SMALL
00340147	Eazy Tray MEDIUM
00340148	Eazy Tray LARGE

# DABLIU MIXING CUPS

Silicone mixing cups, indicated to mix resins. Five different sizes available: 5-7 -10-25-75 ml



<u>00340270</u>	Silicone Mixing Cups 5 ML
<u>00340271</u>	Silicone Mixing Cups 7 ML
<u>00340272</u>	Silicone Mixing Cups 10 ML
<u>00340273</u>	Silicone Mixing Cups 25 ML
<u>00340274</u>	Silicone Mixing Cups 75 ML
<u>00340275</u>	Silicone Mixing Cups SET 5pcs 5-7-10-25-75 ML

# DABLIU CARPET

BACTERIOSTATIC CARPET

Decontaminating bacteriostatic mat, ideal to obtain a high aspesis degree and cleaning of the dental clinic environment. Each mat consists in 30 films and each sheet measurer 115 x 60 cm. Each package contain 5 mats.



**Package:**  
5 carpet 115 x 60 cm

00100131 | Dabliu Carpet (5x30 sheets)

# DABLIU PROTECTOR

BARRIER FILM

Pre-cutted, disposable adhesive film, 45 micron thick, stripping length 15 cm/ 1200 sheets roll in each dispenser. Color: transparent and blue.

**Package:**  
12 rolls /box



00100120 | Dabliu Protector Blue  
00100121 | Dabliu Protector Clear

# ANTISTATIC HEALTH CLOGS

THE PLEASURE OF LIGHTNESS

Our clogs are made from E.V.A. (Ethylene, Vinyl, Acetate), which is an injected and expanded, cross-linked plastic material. This material offers good mechanical and physical property together with a nice tactile and aesthetic aspect.



## Some characteristics:

- Resistance to deterioration in time
- Resistance to atmospheric agents
- Resistance to many oils and solvents
- Highly resistance to use, thanks to elasticity and resistance to ripping and tearing.
- Effect of vibration under shocks and collisions
- Resistance to bacteria, thanks to the cross-linked structure
- Non-toxic
- 100% made in Italy
- Soft, light and comfortable
- Antistatic device
- Machine washable max 60°C
- CE marked



## Available in the following sizes:

n° 36/37 – 38- 39 – 40- 41- 42 -43 -44

## Legend:

- A: light blue
- B: white
- F: fuchsia
- V: green

\*\*\* add colour and sizes code

CODE: 000000\*\*\*

Ex: "Antistatic Health Clogs 36 Light Blue": 00000036A

# SLEEVES

## X-RAY BARRIER ENVELOPES

- Quick tear-off opening
- Soft vinyl for a greater tactile feel
- Fit phosphor plates



00480025 | X- RAY barrier envelopes 188-0 25x38 100 pcs

00480026 | X- RAY barrier envelopes 188-1 44x26 100 pcs

00480027 | X- RAY barrier envelopes 188-2 44x31 100 pcs

00480028 | X- RAY barrier envelopes 188-3 58x23 50 pcs

00480029 | X- RAY barrier envelopes 188-4 82x58 50 pcs

# SLEEVES

## X-RAY

- Fully x-ray unit coverage
- Disposable yet economical
- Soft, pliable polypropylene sleeves
- Hold their position until removed
- Easy application and removal

**Size:** cm 36 large x 23 depth x 76 length



00480012 | X-Ray Sleeves 250 pcs

# SLEEVES

## SUCTION TUBE

- Pre-cut opening for easier application
- Made from very supple poly plastic for a greater tactile feel.

**Size:** 36 x 6,5



00480013 | Suction tube sleeves 500 pcs

# SLEEVES

## CHAIR & HEADREST

### Chair Sleeves

- Eliminate the need for separate headrest covers;
- Avoids the need for harmful disinfectants and aerosols prolonging the life of dental units and stools

**Size:** 200 cm x 72 cm / 36 cm x 25,5 cm

### Headrest Sleeves

- Polypropylene sleeves
- Reduces possible cross-contamination risks

**Size:** 36 cm x 25,5 cm



00480022 | Chair Sleeves 125 pcs



00480019 | Headrest Sleeves 500 pcs



# SLEEVES

## X-RAY SENSOR AND MOUTH MIRROR

Made from very supple poly plastic for greater tactile feel.

Size:

220 mm x 38 mm



00480011 | X-Ray sensor and mouth mirror sleeves 500 pcs

# SLEEVES

## X-RAY SENSOR

Made from very supple poly plastic for greater tactile feel.

Size:

38/58 mm x 220 mm



00480010 | X-Ray sensor sleeves 500 pcs.

# SLEEVES

## INTRAORAL CAMERA

Very supple poly plastic for greater tactile feel.  
Easy to apply

Size: □

23 mm x 120 mm



00480016 | Intraoral Camera Sleeves 500 pcs

# SLEEVES

## T-STYLE LIGHT HANDLE

- Universal sizing to fit most T-style handles
- Full handle coverage
- Disposable yet economical
- Soft polypropylene
- Easy application and removal □

Size:

Universal



00480014 | T-Style light handle 500 pcs

# SLEEVES

## PEN TYPE CURING LIGHT

Made from very supple poly plastic  
Easy to apply

### Size:

320 mm x 84 mm x 26 mm



00480015 | Pen type curing light sleeves 500 pcs

# SLEEVES

## CURING LIGHT

Cover curing light from tip to handle

### Size:

80 mm x 270 mm x 208 mm



00480017 | Curing light sleeves 500 pcs

# SLEEVES

## SYRINGES

- Pre-cut opening for easier application
- Ideal for 3-way syringes, saliva ejectors etc.

### Size:

26 cm x 6,5 cm



00480018 | Syringe sleeves 500 pcs

# SLEEVES

## HANDPIECE

- Disposable sleeves for most handpieces, scalers, made from soft pliable poly plastic.
- Easy to apply.



00480020 | Low speed sleeves 3,75x23 500pz.

00480021 | High speed sleeves 2,5x23 500pz.

# DABLIU LAB

POLYCONDENSATION VINYPOLYSILOXANE

Technical professional use

**Usage:** indicated for laboratory works, as support for provisory prosthesis.



DABLIU LAB 75



DABLIU LAB 85



DABLIU LAB 95

## Characteristics:

- Better elastic properties
- High fluency index
- Less distortion
- Heat resistant (110°C)
- Can be activated with both gel and liquid catalyst
- No greasy surface

0000314	Dabliu Lab 75 5kg
0000315	Dabliu Lab 85 5kg
0000316	Dabliu Lab 95 5kg
0000312	Dabliu Catalyst Paste
0000313	Dabliu Catalyst Liquid

## Characteristics

Iso 4823 - ADA 19			
	Dabliu Lab 75	Dabliu Lab 85	Dabliu Lab 95
Mixing ratio	100:1	100:1	100:1
Hardness after 24H	75 Shore	85 Shore	95 Shore
Linear shrinkage	< 1,0	< 1,0	< 1,0
Mixing time	1 min.	1 min.	20-30 sec.
Total setting time	5 min.	5 min.	4 min.
Color	Neutral	Green	Orange

# DABLIU LAB A+B

MATERIAL FOR HIGH PRECISION IMPRESSION

High viscosity addition silicone

**Usage:** High viscosity addition silicone



DABLIU LAB SUM A+B 70



DABLIU LAB SUM A+B 85



DABLIU LAB SUM A+B 90

## Characteristics:

- High fluency index
- Heat resistant (110°C)
- Optimal elastic recovery
- Less distortion
- Dimensional stability
- Doesn't crumble under drilling
- No greasy surface

0000327	Dabliu Lab 70 5+5kg
0000328	Dabliu Lab 85 5+5kg
0000329	Dabliu Lab 90 5+5kg

## Characteristics

Iso 4823 - ADA 19			
	LAB A+B	LAB A+B	LAB A+B
Hardness	70 Shore	85 Shore	90 Shore
Mixing time	60 sec.	60 sec.	60 sec.
Working time	120 sec.	120 sec.	120 sec.
Setting time	180 sec.	180 sec.	180 sec.
Mixing ratio	1:1	1:1	1:1
Strain in compression	0,2 %	0,2 %	0,2 %
Elastic recovery	99,4 %	99,4 %	99,4 %

# FLUENT CLEAR FLUID

IMPRESSION REALIZATION FOR MASK



Transparent silicon Automix **FLUID**  
for professional use in dental clinic.

**Usage:** Indicated in the realization masks for the moulding of light-curing materials in technique with use of transparent flask thanks to a very fluid and flowing material.

**Indication FLUID:**

- flowability calibrated for use in flask
- high index of fluency
- very high transparency level
- compact and smooth surfaces
- calibrated tear strenght
- high final hardness (70 Sh.A)
- dimensionally stable in time

0000318 | Fluent Clear - White 2 cartridges 50gr.

# FLUENT GUM & SCAN

IMPRESSION GINGIVAL REPRODUCTION



Addition vinylpolysiloxane for gingival reproduction.

**Usage:** Indicated for implant prosthetics for direct method realization of gingival reproduction in pink vinylpolysiloxane easy to trim, for the determination of the shape of the cervical contours and the emergence profile in fixed or mobile prosthesis

**Characteristics:**

- trimming after 10 minutes from extrusion
- high index of flowability
- opaque colour which give it a natural appearance like real gingiva
- elasticity comparable to natural tissue
- balanced tear strenght
- high dimensional stability
- no greasy surface

## Characteristics

Hardness	75 Shore
Iso 4823 - ADA 19	
Mixing time	Automix
Working time	60 sec.
Time molding	6 min.
Time under drill	10 min.
Strain in compression	0,1 %
Elastic recovery	99,8 %
Reproduction of detail	2 micron

0000320 | Fluent Gum & Scan



## DUPLEX LAB 14 - 25

Technical professional use

Usage: used for technical duplication



DUPLEX LAB 14



DUPLEX LAB 25

0000321 Duplex Lab 14 1+1kg  
0000323 Duplex Lab 25 1+1kg

### Characteristics:

- Fast light curing
- High fluency index
- Heat resistant (110°C)
- Optimal elastic recovery
- Precise reproduction of the finest details
- Dimensional stability
- Resistance to splitting
- Compatible with all materials (gypsum, resins)

### Characteristics

Iso 4823 - ADA 19

	Duplex Lab 14	Duplex Lab 25
Mixing ratio	1:1	1:1
Hardness after 24H	14 Shore	25 Shore
Linear shrinkage	< 0,05%	< 0,05%
Mixing time	5 min.	5 min.
Polymerization time	25 min.	25 min.
Total setting time (23°C)	30 min.	30 min.
Reproduction of detail	2 micron	2 micron
Color	Pink	Orange

## DUPLEX ANGLE 12 - 22 - 27

Technical professional use

Usage: laboratory model duplication



DUPLEX ANGLE 12



DUPLEX ANGLE 22



DUPLEX ANGLE 27

### Characteristics:

- Fast light curing
- High fluency index
- Heat resistant (110°C)
- Optimal elastic recovery
- Precise reproduction of the finest details
- Dimensional stability
- Resistance to splitting
- Compatible with all materials (gypsum, resins)

0000324 Duplex Angle 12 1+1kg  
0000325 Duplex Angle 22 1+1kg  
0000326 Duplex Angle 27 1+1kg

### Characteristics

Iso 4823 - ADA 19

	DUPLEX 12	DUPLEX 22	DUPLEX 27
Mixing ratio	1:1	1:1	1:1
Hardness after 24H	12 Shore	22 Shore	27 Shore
Linear shrinkage	< 0,05%	< 0,05%	< 0,05%
Mixing time	5 min.	5 min.	5 min.
Polymerization time	20 min.	20 min.	20 min.
Total setting time (23°C)	25 min.	25 min.	25 min.
Reproduction of detail	2 micron	2 micron	2 micron
Color	Blue	Light Blue	Green

# DABLIU WAX

Multitasking red wax in sheets, for bases, fixing and waxing in removable prosthesis. Thickness 1,5 mm.  
Box 500g.

00170006 | Dabliu Wax 500 gr.



# DABLIU PUMICE

Natural Pumice for the final polishing of prosthesis works.  
5kg plastic drum.

00280030 | Natural powder pumice 5 kg



Red Top



# DABLIU STONE IV - EH

Extra Hard Gypsum type IV for high precision models and dies. Added with synthetic resin, thanks to its high dimensional stability, resistance to abrasion and plasticity, there is no risk of chipping during its working. It is compatible with all impression materials. 5 kg plastic drum.

## Technical characteristics:

Mixing ratio 100g : 20 ml  
Mixing time M/A: 60 sec. / 30 min.  
Working time: 8 min.  
Final setting time: 12 min.  
Removal time: 30 min.  
Setting Expansion: <0,10% (2h)  
Compressive strength: >50 MPa (1h)

00280018 | Extra hard gypsum Type IV Golden brown 5 kg



Brown Top



# DABLIU STONE MODEL IV

Extra fine Hard Gypsum Type IV, with high precision for models.

## Technical characteristics:

Mixing ratio 100g : 22-23 ml  
Mixing time M/A: 60 sec. / 30 min,  
Working time: 6 min.  
Final setting time: 8 min.  
Removal time: 30 min.  
Setting Expansion: <0,15% (2h)  
Compressive strength: >45 MPa (1h)

00280010 | Extra Hard Gypsum Type IV 5kg

00280020 | Extra Hard Gypsum Type IV 25 Kg



Yellow Top



## DABLIU STONE ORTHO III

A hard type III gypsum, for dentist's surgery and archive orthodontic models. Due to its good fluidity it can be used with acrylic resins at hot and cold temperatures. Compatible with all impression materials. 15 kg plastic drum.

### Technical characteristics:

Mixing ratio 100g : 26-28 ml  
 Mixing time M/A: 60 sec. / 30 min.  
 Working time: 6 min.  
 Final setting time: 8 min.  
 Removal time: 30 min.  
 Setting Expansion: <0,20% (2h)  
 Compressive strength: >30 MPa (1h)



White Top



00280026 | Orthodontic white gypsum – Type III 15kg

## DABLIU STONE ORTHO IV

An extra hard type IV gypsum used for precision models in fixed orthodontics. It is compatible with all impression materials. 5 kg plastic drum.

### Technical characteristics:

Mixing ratio 100g: 22-23 ml  
 Mixing time M/A: 60 sec. / 30 min.  
 Working time: 6 min.  
 Final setting time: 8 min.  
 Removal time: 30 min.  
 Setting Expansion: <0,10% (2h)  
 Compressive strength: >45 MPa (1h)



White Top



00280017 | Extra Hard White Gypsum - Type IV 5 Kg

## DABLIU STONE MODEL III

Hard gypsum, type III, for high precision for antagonist models and flack. Extra fine-grained it assures smooth surfaces. The model results dense and inalterable during use. It is compatible with all impression materials.

### Technical characteristics:

Mixing ratio 100g: 28-30 m  
 Mixing time M/A: 60 sec. / 30 min.  
 Working time: 6 min.  
 Final setting time: 8 min.  
 Removal time: 30 min.  
 Setting Expansion: <0,20% (2h)  
 Compressive strength: >30 MPa (1h)



Yellow Top



00280041 | Yellow gypsum Type III 15 kg

00280021 | Yellow gypsum Type III 5 kg

## DABLIU STONE MODEL II

A soft type II gypsum, for model bases, for occlusal moulds and laboratory models. Extremely fine, free of impurities. 25 kg paper bag.

### Technical characteristics:

Mixing ratio 100g: 55-60 ml  
 Mixing time M/A: 60 sec. / 30 min.  
 Working time: 5 min.  
 Final setting time: 8 min.  
 Removal time: -  
 Setting Expansion: <0,30% (2h)  
 Compressive strength: >10 MPa (1h)



White Top



00280048 | Soft gypsum Type II 25kg

# Mimetic

BRACKETS

Sapphire brackets, unique and completely transparent, that totally blend with the dental structure.



# Eclipse

ARCH



Ni-Ti archwires created to remain active until its full form is achieved. Their surface allows easy wire sliding

# Ortho Control



S&TCR are two fixed orthopaedic devices functioning with proprioceptive stimulus, that thanks to their 24 hours-a-day action, they get the therapeutic result, independently from the patient cooperation.

Ortho Control are preformed appliances for intercepting orthodontics, indicated for growing patients, to prevent malocclusions.



**S&TCR**  
Second & Third  
Class Resolver

*Perfect*



*LINE*

*Made in Italy*



Complete line of orthodontic instruments, resistant and durable. Cutters are TC edge endowed.